

OCULUS | Mesotest II



INSTRUCTION MANUAL

Examination of mesopic vision and glare sensitivity

Notes on this Instruction Manual

To ensure safe operation, it is essential that you use the device correctly. For this reason you should familiarize yourself thoroughly with the contents of this instruction manual before operating the device. In particular, pay attention to the safety instructions.

The Mesotest II is used for testing the mesopic visual acuity of the human eye.

Due to ongoing development, the diagrams shown may depict minor changes to the actual device delivered.

If you have any queries or would like additional information about your device, do not hesitate to call or send us a mail or a fax. Our service team will gladly assist.



OCULUS is certified according to DIN EN ISO 13485, setting high standards of quality for the development, manufacture, quality assurance and service of the entire range of products.

Table of Contents

1	Scope of Delivery.....	1
2	Graphic Symbols on the Equipment	2
3	Information about the Documentation	2
3.1	Structure of the Documentation.....	2
3.2	Software Version	3
4	Safety Instructions.....	3
4.1	About this Manual	3
4.1.1	Used Graphic Symbols.....	4
4.2	Safety Instructions for Use	4
5	Intended Use.....	8
6	Device Description	9
6.1	Mesotest II Functionality	11
6.2	Twilight (Mesopic) Vision	14
7	Initial Operation	16
7.1	Before Initial Operation	16
7.2	Information About the Room in which the Device is Set Up	16
7.2.1	Insert Cover	17
7.3	Set up and Connection	17
7.4	Software Installation.....	19
7.5	Switching the Device On and Off	19
8	Before Conducting the Examination	20
8.1	Prepare the Unit for Routine Operation.....	20
8.2	Prepare the Patient for the Examination.....	20
8.2.1	Patient Adaptation	21
9	Examination Using the Control Pad.....	22
9.1	Control Pad.....	22
9.2	Filling In the Test Sheet.....	25
9.2.1	Starting the Examination	25
9.3	Test Mesopic Vision	26
9.4	Test Mesopic Vision With Glare	26
9.4.1	Test for Night Myopia (optional).....	27
9.4.2	Monocular Examinations (optional).....	27
10	Conducting the Examination with a Computer/Laptop	28
10.1	Using the Examination Menu	28
10.1.1	View the Patient's Information and Enter any Visual Aid.....	29
10.1.2	Marking the Results in the Test Image Display.....	30
10.1.3	Present the Next Optotype	30
10.1.4	Swing in Minus Lenses.....	31
10.2	Change Settings.....	31
10.3	Conduct the Examination	36

10.3.1	Starting Patient Data Management	37
10.3.2	Entering a New Patient.....	38
10.3.3	Selecting an Existing Patient	38
10.4	Select and Start Exam.....	39
10.4.1	Conduct the Programm "Mesopic vision bin".....	40
10.4.2	Conduct the Programm "Mesopic vision mono".....	42
10.4.3	Test for Night Myopia (optional).....	44
10.5	Display and Print Results.....	45
10.6	Finishing the Exam.....	46
10.7	Load Existing Examination.....	46
11	Editing the Examination Programmes.....	47
11.1	Compilation of a Custom Selection List of Examinations	47
11.2	Create A New Examination Programme	48
11.3	Edit Examination Programme	51
12	Patient Data Management.....	52
12.1	Rename Patient Data	52
12.2	Exporting Patient Data	53
12.3	Importing Patient Data.....	54
12.4	Data Backup.....	55
12.4.1	Backup data.....	56
12.4.2	Reconstruct data.....	56
12.4.3	Automatic backup	57
12.5	Change settings	57
12.5.1	"Main" tab	58
12.5.2	"Devices" tab.....	60
12.5.3	"Import/Export" tab.....	62
12.5.4	"Email" tab.....	64
12.5.5	"Interface" tab	65
12.5.6	"Smartcardreader" tab.....	66
12.5.7	"Misc" tab	67
12.5.8	"Hecht" and "DICOM" tabs.....	68
13	Cleaning, Disinfection and Maintenance	68
13.1	Unplug the Equipment.....	68
13.2	Cleaning	69
13.3	Disinfection	69
13.4	Care and Maintenance.....	70
13.4.1	Exchange Viewer Enclosure	71
13.4.2	Replacing the Fuse	72
13.4.3	Change the Lamp at the Panel Illumination	73
14	Troubleshooting.....	74
15	Transport and Storage	75
15.1	Transport and Storage Information.....	75
15.2	Disassembly and Packing	76

16 Disposal.....	76
17 Terms of Warranty and Service.....	77
17.1 Liability for proper function or damages.....	77
17.2 Address of the manufacturer and service department.....	78
18 Declaration of Conformity.....	79
19 Technical Data	80
20 Appendix.....	82
20.1 Electromagnetic Compatibility.....	82
20.2 Guidance and Manufacturer's Declaration - Electromagnetic Emmissions for the Mesotest II.....	83


1 Scope of Delivery

Components	Order Number
■ Mesotest II	62800
or	
■ Mesotest II konkav (minus lenses can be swung in), with eye occluder	62801
■ Mains cable	05200320
■ Test block (100 sheets)	62850
■ Dust protection cover	026010005001
■ Cover for back pane of the viewer	016280000002
■ Instruction Manual Mesotest II	G/62800/xxxx/en
Optional:	
■ Control pad with LCD module and panel illumination	62811
■ Software module for Windows®-compatible computers for patient data acquisition and for control of the Mesotest II, incl. interface cable, USB 2.0 adapter and floating license key	62803
■ Carrying bag	62802

We reserve the right to change the scope of delivery in line with ongoing technical development.

- ➔ If you find transport damage upon delivery, immediately file a claim with the transport company.
- ➔ Have the damage confirmed on the bill of lading so that an orderly handling of the complaint for damages is possible.
- ➔ Keep the packaging in a safe place. You can then ship or transport the unit in the proper manner for any servicing or repairs that may arise. You can thus avoid unnecessary damage and costs.

2 Graphic Symbols on the Equipment

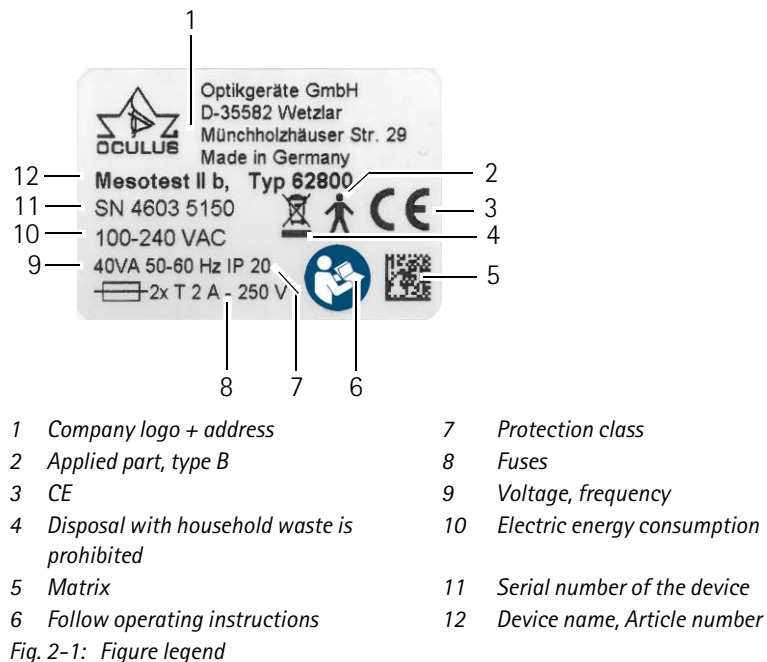

3 Information about the Documentation

3.1 Structure of the Documentation

A folder containing a set of documentation is supplied with your Mesotest II:

- **Instruction Manual:** The design of the unit is described in detail in this document. The instruction manual also gives you general information about working with the Patient Data Management system and all safety-related instructions for use of the Mesotest II.
- **Software Installation:** The introduction to the Software Installation describes how to install the Mesotest II software and the associated drivers.

If you work with a floating license key, the respective instruction manual explains how you can use the Mesotest II within a network.

3.2 Software Version

The instruction manual describes the following versions of the Mesotest II software as well as the Patient Data Management system:

Laptop/Computer:

- Mesotest II software: Version 1.0.4.x
- Patient Data Management Version 6.08



Note:

- The Patient Data Management software version is displayed on the „Settings“ screen within the Patient Data Management ([sec. 12.5, page 57](#)).
 - The software version of the Mesotest II programme is shown on the “Change Settings” screen within the Mesotest II programme ([fig. 10-11, page 35](#)).
-

4 Safety Instructions

4.1 About this Manual

- ➔ Carefully read through the Instruction Manual.
- ➔ Keep the Instruction Manual in good condition near the device.
- ➔ If the scope of delivery includes a netbook: Read the operating manual for the netbook.
- ➔ Observe the legal regulations with regard to accident prevention.

4.1.1 Used Graphic Symbols



Attention

Denotes a potentially hazardous situation which can easily result in minor physical injury or property damage.



Note:

Identifies situations which can result in property damage, or denotes application information and useful or other important information.



Identifies important information about the product and its use, which require special attention.

- > This symbol denotes menu paths and screen shots. Example for starting an examination:
Mesotest II > Examination > Load
which means:
 - ➔ Select the "Examination" menu from the menu bar.
 - ➔ Select the menu item "Load".

4.2 Safety Instructions for Use



Attention

Personal injury or property damage due to improper operation

- ➔ Observe the following safety instructions.

Personal injury or property damage due to equipment modifications that could jeopardize safety

- ➔ This equipment may not be modified without the permission of the manufacturer.

Instructions for operating personnel

- ➔ Make certain that the Mesotest II is used exclusively by personnel that have the training and practical experience to safely and properly operate the equipment.

Transport and storage instructions

Refer to the notes in *sec. 15, page 75*.

Instructions for setup and connection

- ➔ Do not use or store the Mesotest II in rooms that are humid.
- ➔ Keep the Mesotest II away from water that may drip, surge or splash on it, and make sure that no liquids can enter the Mesotest II. Do not place any containers holding liquids in the vicinity of the Mesotest II.
- ➔ Germany: Only operate the Mesotest II in rooms used for medical purposes after they have been set up according to the VDE Regulations 0100-710.
- ➔ Do not operate the equipment included in the packing list in explosive environments, in the presence of combustible narcotic agents or volatile solvents such as alcohol, benzene, etc.
- ➔ Make sure that the electrical plug of the Mesotest II is easily accessible so that it can be unplugged easily.
- ➔ Do not use excessive force when connecting the electrical plug. If a connection is not possible, check whether the plug fits the jack. If you find damage to the plug connector, have the damage corrected by our service department.

Information about patient environment

Patient environment is the area where patients can come into contact with any part of a medical electrical equipment (ME equipment) or with another person being in contact with the ME equipment.

Use devices in the patient environment which are in conformity with the IEC 60601-1 standard. If you use a multiple socket extension cord or a device which do not comply with the IEC 60601-1 standard, use an isolation transformer

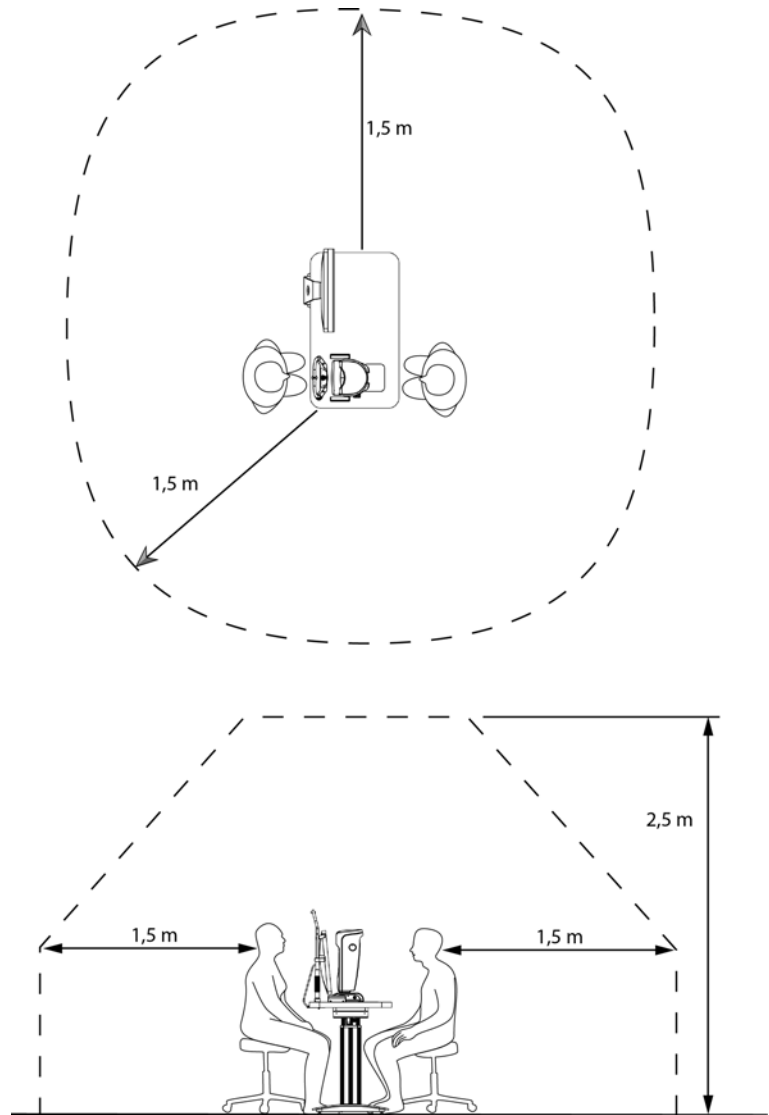


Abb. 4-1: Patient environment

Instructions to use a ME system

The Mesotest II and a connected computer are a medical electrical system (ME system) according to DIN EN 60601-1. If you want to connect further devices, for example a printer, these devices will be part of the ME system.

- ➔ Make sure that all devices of the ME system are in accordance with IEC 60601-1 or IEC 60950-1.

Instructions for operation

- ➔ Never put a damaged Mesotest II into operation.
- ➔ Only operate the Mesotest II using original accessory parts supplied by us, and when the device is in technically perfect working order.
- ➔ Before initial operation: Let OCULUS or an authorized dealer train you in the operation of the Mesotest II.
- ➔ Only operate the equipment after you have read and understood the Instruction Manual.
- ➔ Do not touch the patient and the device at the same time.
- ➔ Make sure that the device cannot tip over by leaning against it or sitting on it.
- ➔ Do not cover the ventilation openings.

Instructions for maintenance

To ensure satisfactory and reliable operation, we recommend that you have the Mesotest II checked every two years by our service department or an authorized dealer. If an error occurs which you are unable to correct, label the Mesotest II as "out of order" and contact our service department or an authorized dealer.

Instructions for disassembly and disposal

- ➔ When disconnecting electrical connections, pull on the respective plug instead of the cable itself.
- ➔ Dispose of the equipment in compliance with the corresponding legal requirements.

Note on electrical safety



Attention

Risk of personal injury or damage to property due to an incorrect level of safety

Connecting the Mesotest II with non-medical electrical equipment (for example data processing equipment) to a medical electrical system must not result in a patient safety level below that prescribed by DIN EN ISO 60601-1. If making this connection leads to the leakage current threshold being exceeded, protective measures including a circuit breaker must be in place.

- ➔ Ensure that connections with non-medical devices are made correctly.
- ➔ Use only a computer that meets the specifications given in this instruction manual, [sec. 19, page 80](#).



Attention

Use of a multiple socket extension cord

Risk of personal injury or material damage caused by unsafe multiple socket extension cord

If you are using a multiple socket extension cord to connect the Mesotest II to the power supply, you must heed the following information:

- ➔ Use an extension cord strip that complies with the requirements of DIN EN ISO 60601-1: 20005, section 16.
- ➔ Do not place the multiple socket extension cord on the floor.
- ➔ Do not use more than one multiple socket extension cord.
- ➔ Plug only the Mesotest II and the computer that is being used with the unit (if applicable) into the multiple socket extension cord.

If you are using a multiple socket extension chord it has to be supplied with a isolation transformer.

If you are using a new computer for the Mesotest II:

- ➔ Check the electrical safety. Call OCULUS Service for this purpose.

Electromagnetic Compatibility (EMC) / cables

Risk of personal injury or damage to property due to electromagnetic interference

Portable and mobile RF communications equipment can affect medical electrical equipment [sec. 20, page 82](#).

- ➔ Make sure that portable and mobile RF communications equipment do not cause interference.
- ➔ Recommendation: Maintain a minimum distance of 4 m. If the distance is shorter, you must ensure that the Mesotest II functions correctly.

5 Intended Use

The Mesotest II is intended for the use described in this Instruction Manual.

The Mesotest II is used for testing the mesopic visual acuity of the human eye.

- ➔ Heed the safety instructions listed above.

Contraindication

None known

6 Device Description



1 Viewer enclosure


2 Device viewer

3 Housing

4 Yoke spring

5 Air vents

Fig. 6-1: Front view of the Mesotest II

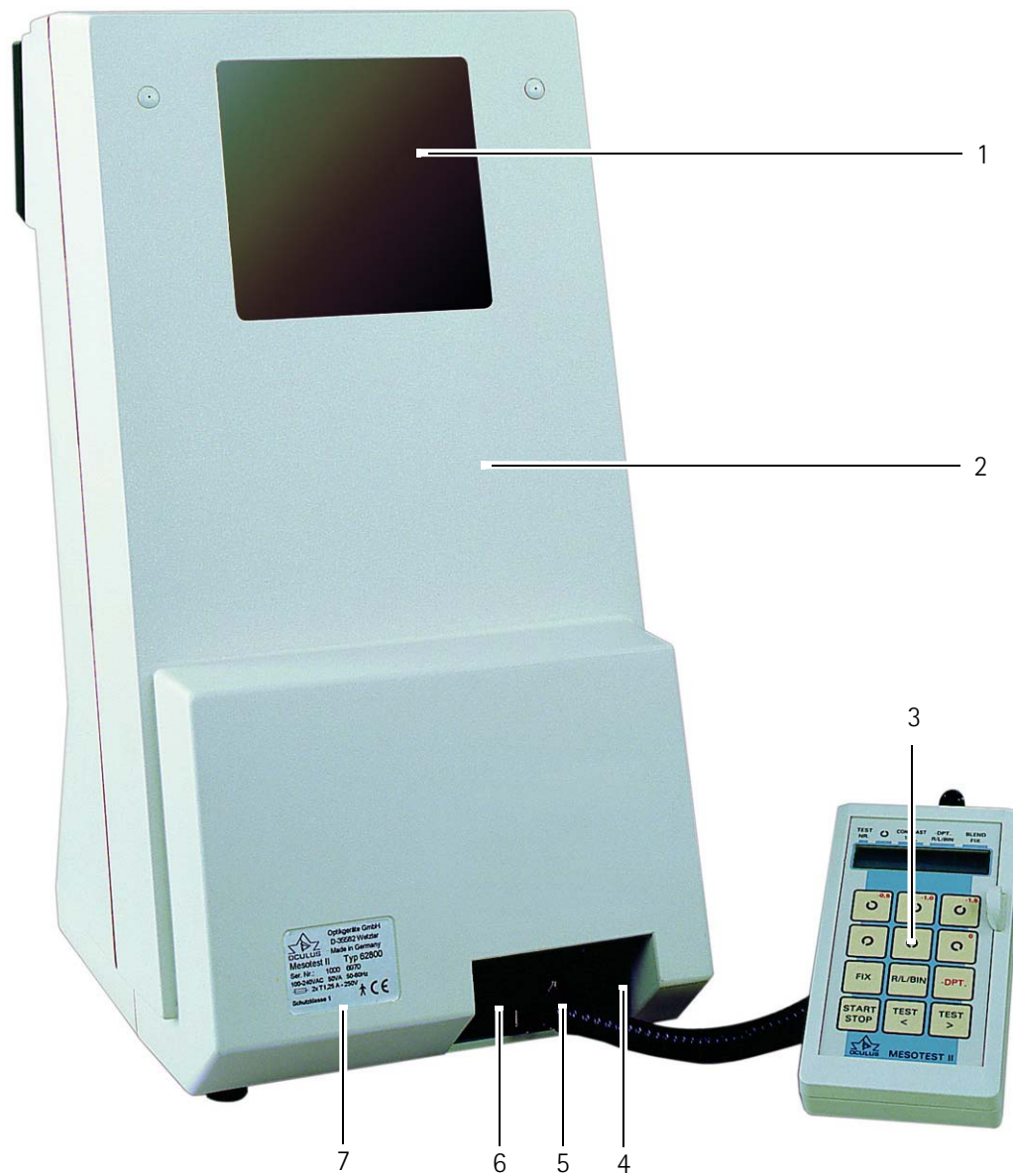

- 1 Back of the viewer with grey pane
- 2 Housing
- 3 Control pad
- 4 Device plug with fuse holder

- 5 Serial interface port
- 6 On/Off switch
- 7 Nameplate

Fig. 6-2: Back view of the Mesotest II

6.1 Mesotest II Functionality

Description

The Mesotest II is a PD-independent, freespace instrument for testing mesopic vision and glare sensitivity.

Freespace means that the tests are conducted under natural visual conditions in an open space.

The device myopia is more or less eliminated, as the eyes accommodate and converge in the same way as for natural vision.

The flexible viewer enclosure facilitates relaxed viewing and provides an optimal fit for every patient. The enclosure completely shuts out all extraneous light, thus allowing correct examinations to be conducted, even in just moderately dark room conditions.

The air vents in the enclosure prevent the patient's glasses or the unit's viewer glass from fogging.

Use

In Germany: The Mesotest II is used for testing mesopic vision and glare contrast sensitivity as required by driver's license insurance regulations and it meets the requirements of DIN 58220, Part 7.

Testing of mesopic vision and glare contrast sensitivity is an important supplement to photopic vision testing, especially for the assessment of night driving ability. Pre- and post-operative testing of the mesopic vision and glare sensitivity of patients with intraocular lenses, or those who have undergone refractive surgery is recommended, as well as for prescription of tinted lenses.

The device is dimensioned such that the patients' differing inter-pupillary distance does not affect the results.

Setup

The focussing screen of the test field is viewed through a lens via a semi-reflector viewer mirror. The optotypes are presented in front of this focussing screen. Their virtual image appears 5 meters away from the eye. The optotypes are vapour-deposited on a glass substrate.

The glare light is located at the left of the test field at a visual angle of 3°. Depending on the test that has been selected, this is activated automatically.

The glare light simulates the glare situation of an oncoming vehicle with dimmed headlights. The glare illuminance at pupil level is 0.35 lux.

White light emitting diodes are used as the light source. It is thus unnecessary to change lamps

The ambient light and the glare light are illuminated separately by one light emitting diode each.

Controls

You can operate the Mesotest II either via the control pad or via a computer/laptop.

Control Pad: The integrated panel illumination allows testing to be done in a darkened room, also see [sec. 9.1, page 22](#).

Operation via computer/laptop: You can also operate the Mesotest II via a computer/laptop. You then need the Mesotest II software (available as an optional extra), also see [sec. 10, page 28](#).

Testing with Landolt rings

Testing is conducted under natural vision conditions, without device accommodation.

Testing is done with a Landolt ring of the visual acuity level 0.1 in acc. with DIN 58220, Part 1 (in Germany). The Landolt ring can be presented in six different positions. Four contrast levels without glare (ambient luminance 0.032 cd/m²) and four contrast levels with glare (ambient luminance 0.10 cd/m²) are tested.

The optotypes can be presented in six different positions, whereby one position is prespecified at the beginning of the test.

Four different contrast levels are available:

1:23 / 1:5 / 1:2.7 / 1:2

The contrast level 1:23 is the strongest contrast, i.e. this level is the easiest to see.

Contrast 1:23 is the ratio of

Ambient luminance: Luminance of the optotype.

Glare light

The glare light is located on the test field.



- 1 Landolt ring
- 2 Ambient light
- 3 Glare light

Fig. 6-3: Glare light

The glare light simulates the glare situation of an oncoming vehicle with dimmed headlights in the twilight (glare sensitivity). The glare illuminance at pupil level is 0.35 lux.

White light emitting diodes are used as the light source. It is thus unnecessary to change lamps.

The ambient light and the glare light are illuminated separately by one light emitting diode each.

A stand-by function switches off the illumination after a certain time if no entries are made by the operator.

Fixation aid

A projection device for a fixation aid is located over the test field. Two red dots can be projected onto the test field with this device, if necessary.

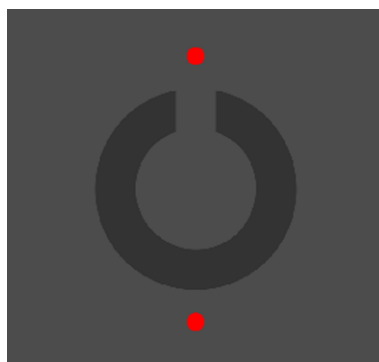


Fig. 6-4: Fixation aid

If the patient looks directly into the glare light when the switch-over to mesopic vision is made, he will no longer be able to see the optotypes, to begin with. This can lead to a delay in the examination procedure, as he needs some time (approx. 10-30 seconds) to readapt.

Test procedure

Control Pad: The contrast levels are each presented with and without glare. As such, a total of eight consecutive tests are conducted.

Test no.	1	2	3	4	5	6	7	8
Contrast	1:23	1:5	1:2.7	1:2	1:23	1:5	1:2.7	1:2
Glare light	off	off	off	off	on	on	on	on
Prespecified position of the Landolt ring	TL	TR	T	B	BR	B	BL	T
Function	Test mesopic vision				Test mesopic vision with glare			

The visual requirement is deemed to be met when the patient has recognised at least 3 of 5 positions of the Landolt ring (60% criterion). Mark this accordingly on the test sheet.

Computer/Laptop:

With the examination program „Mesopic vision bin“, you can examine mesopic vision binocularly.

It consists of 10 test steps.

Test 1 and Test 6: Facilitate the fixation, as here the fixation marks are displayed

Test 2-5: Are used to test mesopic vision

Test 7-10: Are used to test mesopic vision with glare.

With the examination program „Mesopic vision mono“ (optional with Mesotest II konkav), you can examine mesopic vision monocularly, for the right and for the left eye. It consists of 20 test steps.

Test 1, 6, 11 and 16: Facilitate the fixation, as here the fixation marks are displayed.

Test 2-5 and 12-15: Are used to test mesopic vision

Test 7-10 and 17-20: Are used to test mesopic vision with glare.

5 Landolt rings are presented for each contrast level. The visual requirement is deemed to be met when the patient has recognised at least 3 of 5 positions of the Landolt ring (60% criterion). Mark this in the software accordingly ([sec. 10.1.2, page 30](#)). The results are saved and you can print them out. You can also compile your own tests.

6.2 Twilight (Mesopic) Vision

Tests for determining mesopic visual acuity and glare contrast sensitivity were developed, among other things, to test the visual functions that a driver must have for driving at night.

How do we see in twilight?

The transition between daytime (photopic) vision and night (scotopic) vision is referred to as twilight vision (mesopic vision).

Twilight vision has considerable differences compared to daytime vision.

Daytime vision and cone vision (also known as photopic vision) are the terms used for luminances greater than 3 cd/m², night vision and rod vision (also known as scotopic vision) for luminances below 0.03 cd/m², and mesopic vision for luminances in the range between 0.03 and 3 cd/m².

Testing mesopic vision	To adapt the test conditions to the luminance conditions prevailing on roads at night, the Landolt ring is presented on a test field with a luminance of 0.032 cd/m², and for the test with glare, on a test field with a luminance of 0.1 cd/m². During the function test, unlike in the visual acuity test, the size of the Landolt ring is not reduced, but rather the contrast that it has to its surroundings. The size of the Landolt ring corresponds with a visual acuity of 0.1. Four different contrast levels are available for reducing the contrast. ■ 1:23 ■ 1:5 ■ 1:2.7 ■ 1:2 The contrast level 1:23 is the strongest contrast, i.e. this level is the easiest to see.
Testing mesopic vision with glare	Testing of mesopic vision can also be conducted with additional glare. The glare comes from the left, just as it would in a road traffic situation.
Thresholds according to DOG (German Ophthalmological Society)	The following thresholds for assessment of fitness for driving in road traffic are recommended by DOG and the German Association of Ophthalmologists: ■ Classes D, D1, DE, D1E: Contrast 1:2,7 ■ Classes C, C1, CE, C1E and taxi drivers: Contrast 1:5 ■ Classes A, A1, A2, B, BE, AM, L and T: Contrast 1:23 A stop criterion similar to that described in DIN 58 220 (visual acuity test) must be established as the test criterion: 3 out of 5 of the different optotypes presented must be correctly recognised. This means that for the critical contrast levels 1:23 to 1:2.7, at least 5 different optotypes must be presented.

7 Initial Operation

7.1 Before Initial Operation

Before the Mesotest II is operated for the first time, you must complete the following actions:

- Set up
- Connection
- If necessary, install the software



Attention

Risk of incorrect measurements/equipment damage due to improper set-up

- ➔ Before first use, make sure the installation and connection of the Mesotest II are completed by our service or by a professional authorized by OCULUS.

7.2 Information About the Room in which the Device is Set Up

The Mesotest II should be set up in a moderately bright or darkened room (room brightness between 2 and 5 lux).



Attention

Interruption of the examination due to overheating

If the unit overheats, it shuts itself off and the test results are lost.

- ➔ Do not cover the ventilation openings.



Note:

Falsified examination results if the room is too bright

If you set up the Mesotest II in a bright room, you must occlude the back pane of the viewer with the supplied cover.

- ➔ Insert the cover, see [sec. 7.2.1, page 17](#)
-

7.2.1 Insert Cover



Occlusion of the back pane of the viewer, however, increases the risk of device myopia.

Take this into account in your examination results.

- ➔ Insert the cover into the slot between the back pane and the housing at one side.
- ➔ Keep hold of the cover at the recessed grips with one hand.
- ➔ With the other hand, bend the cover and let it snap into the opposite slot.



Fig. 7-1: Insert cover

- ➔ Removing the cover :
Pull out the cover at the recessed grips.

7.3 Set up and Connection



Note:

Risk of damage to the device if not handled correctly

- ➔ Do not expose the Mesotest II to any vibrations, shocks, contaminants, moisture, or high temperatures.
- ➔ Handle the optical device with care.
- ➔ Use a power socket that has a flawless earth conductor connection.

The ambient conditions for operation are given in "[Ambient conditions](#)" [auf Seite 80](#).

- ➔ Remove the Mesotest II from the packaging.
- ➔ Place the Mesotest II on a level surface.
- ➔ Set up the Mesotest II so that the power plug is easy to access. That way, you can easily disconnect it from the power supply for any repairs or maintenance work.
- ➔ Plug the control pad into the jack ([fig. 7-2, page 18, item 2](#)).



Attention

Electrical safety hazard


- ➔ Do not use the Mesotest II immediately adjacent to other devices and do not stack it with other devices.
- ➔ If you use a multiple socket extension cord to connect the Mesotest II: Use a multiple socket extension cord that complies with the requirements of DIN EN ISO 60601-1.
- ➔ Do not place the multiple socket extension cord on the floor.
- ➔ Do not use more than one multiple socket extension cord.
- ➔ Plug only the Mesotest II and the computer that is being used with the unit (if applicable) into the multiple socket extension cord.
- ➔ Use a socket with a protective earth connection which is fully operating.



If you are going to operate the Mesotest II via a computer/laptop, connect the latter in the same way as the control pad, [fig. 7-2, page 18, item 2](#).

- ➔ Connect the unit to the mains with the supplied power cable, [\(fig. 7-2, page 18, item 1\)](#).

Make sure that the mains voltage is the same as the voltage specified on the rating plate.

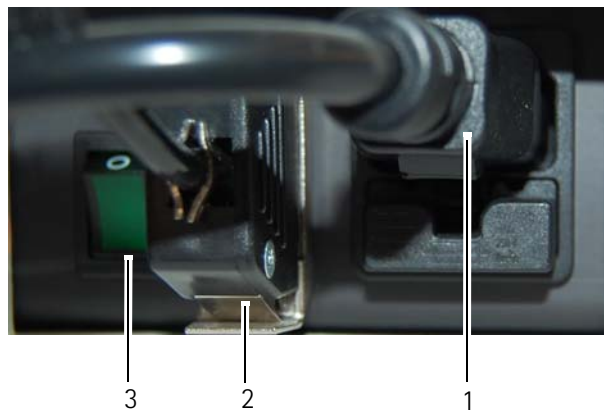

- 1 Mains connection
- 2 Jack for connection of the control pad/laptop/computer
- 3 On/Off switch

Fig. 7-2: Mesotest II Connection

- ➔ Switch on the Mesotest II at the On/Off switch (3).
After a short, automatic self-test, the device is ready for operation. The On/Off button lamp (3) lights up for control purposes. If you are working with the control pad, the LCD module and the panel illumination also light up here.

7.4 Software Installation

Before first use, you must install the Mesotest II software on your laptop or PC. Proceed as described in the [Software Installation](#).

For installation on a computer/laptop, we recommend:

- 2 GB RAM, Windows® 8, free USB port



As of April 2014, Microsoft® no longer provides support for the operating system Windows™ XP. This means that Microsoft® no longer provides any operating and security updates for Windows™ XP.

If you continue to use this operating system after support has ceased, your computer and your data are susceptible to greater security risks. We recommend that you:

- Disconnect your computer from the internet and use it as a stand-alone solution, or
 - Upgrade to the latest Windows™ version
-



You can only conduct the examination with the computer when you use the supplied floating license key and the current version of the device software has been installed on the computer.

- ➔ Plug the floating license key into one of the computer's USB ports. The software then installs itself.
-

7.5 Switching the Device On and Off

- ➔ Switch the Mesotest II on or off at the On/Off switch at the unit ([fig. 7-2, page 18, item 3](#)).

If the Mesotest II is connected with a control pad, the unit starts automatically.

If the Mesotest II is connected with a netbook, laptop or PC, you must start the software, [sec. 10.3.1, page 37](#).

8 Before Conducting the Examination

8.1 Prepare the Unit for Routine Operation

To start the Mesotest II for daily operation, you must:

- Check the connections, [sec. 7.3, page 17](#)
- Check whether the viewer is dirty
- ➔ If impurities are found: Clean the viewer, also refer to [sec. 13.2, page 69](#).
- ➔ Disinfect the viewer enclosure after every examination and the housing whenever necessary, [sec. 13.3, page 69](#).

8.2 Prepare the Patient for the Examination

Prepare the patient for the examination, as follows:

- ➔ Check to be sure that any adaptation disorder has been resolved, [sec. 8.2.1, page 21](#).
- ➔ Check the patient's glasses for damage or dirt.
- ➔ Explain the procedure and the optotypes to the patient. You can show the patient the following figure:

A ring with a gap in it is used as the test symbol (optotype) The gap can appear in eight different positions. Tell the patient where the gap is.	
	Top (12 o'clock)
	Top right (1 to 2 o'clock)
	Right (3 o'clock)
	Bottom right (4 to 5 o'clock)
	Bottom (6 o'clock)
	Bottom left (7 to 8 o'clock)
	Left (9 o'clock)
	Top left (10 to 11 o'clock)

Fig. 8-1: Example: Landolt rings explained



- ➔ Make sure that the viewer enclosure remains in contact with the patient's head during the entire examination. Interfering incidence of light and impairment of the adaptation state associated with it can thus be avoided.

8.2.1 Patient Adaptation



Before the examination, the patient's eyes need about 5 minutes to adapt to the darkness.

When the Mesotest II is set up in a darkened room, this is not a problem. If the unit is set up in a bright room, adaptation of the patient must be done by having him look into the unit for the required length of time prior to starting the examination.

In a darkened room

- ➔ Let the patient's eyes adapt for about 5 minutes.

In a bright room


- ➔ Let the patient look into the unit for about 5 minutes.

9 Examination Using the Control Pad

This chapter explains:

- The control pad, [sec. 9.1, page 22](#)
- How to fill out a test sheet, [sec. 9.2, page 25](#)
- How to start an examination, [sec. 9.2.1, page 25](#)
- How to test mesopic vision, [sec. 9.3, page 26](#)
- How to test mesopic vision with glare, [sec. 9.4, page 26](#)

9.1 Control Pad

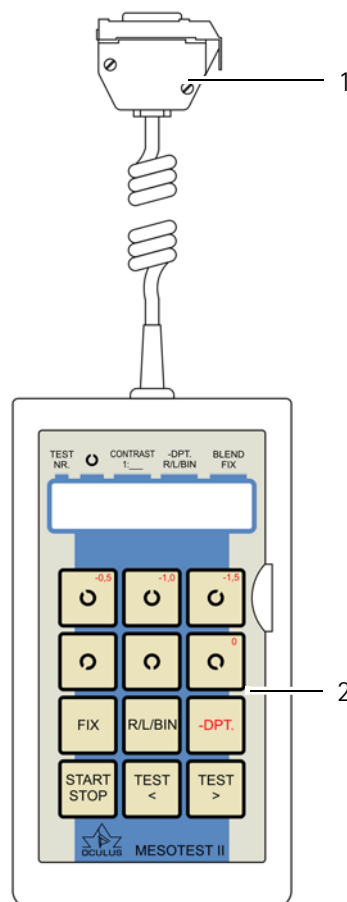

1 Mains connection

2 Control pad

Fig. 9-1: Control pad

Control pad user interface

Button	Meaning
	Landolt ring position: Top Left (LO) 2nd function: interpose -0.5 D minus lenses (option)
	Landolt ring position: Top (O) 2nd function: interpose -1.0 D minus lenses (option)
	Landolt ring position: Top Right (RO) 2nd function: interpose -1.5 D minus lenses (option)
	Landolt ring position: Bottom Left (LU)
	Landolt ring position: Bottom (U)
	Landolt ring position: Bottom Right (RU) 2nd function: reset minus lenses to 0 D (option)
	Turns Fixation Assistance on and off
	Select an eye for examination (option)
	Activates 2nd function (minus lenses) (option)
	Next patient, or stop current tests and return to beginning
	Return to previous test
	Proceed to next test

Control pad user interface

Display	Meaning	
TEST NR. 1 CONTRAST 1: 23 -DPT. R/L/EIN B BLEND FIX 1 LO 23 B	Test no. 1	Display after switch-on
	Optotype: LO	Landolt ring gap at top left
	CONTRAST 1: 23	Presented contrast level 1:23
	-DPT. R/L/BIN: B	Binocular examination, without minus lens
	BLEND FIX: (empty)	Glare and fixation aids are not switched on
TEST NR. 5 CONTRAST 1: 23 -DPT. R/L/EIN B BLEND FIX FB 5 RU 23 B FB	Test no. 5	Fifth programme step
	Optotype: RU	Landolt ring gap at bottom right
	CONTRAST 1: 23	Presented contrast level 1:23
	-DPT. R/L/BIN: B	Binocular examination, without minus lens
	BLEND FIX: FB	Glare and fixation aids are switched on
TEST NR. 7 CONTRAST 1: 2.7 -DPT. R/L/EIN B BLEND FIX B 7 O 2.7 B B	Test no. 7	Seventh programme step
	Optotype: O	Landolt ring gap at top
	CONTRAST 1: 2.7	Presented contrast level 1:2.7
	-DPT. R/L/BIN: B	Binocular examination, without minus lens
	BLEND FIX: B	Glare is switched on and fixation aids are switched off
TEST NR. 6 CONTRAST 1: 5 -DPT. R/L/EIN B 10 B BLEND FIX 6 U 5 B 10 B	Test no. 6	Sixth programme step
	Optotype: U	Landolt ring gap at bottom
	CONTRAST 1: 5	Presented contrast level 1:5
	-DPT. R/L/BIN: B 10	Binocular examination, minus lens -1.0 dpt is swung in (option)
	BLEND FIX: B	Glare is switched on and fixation aids are switched off
TEST NR. 4 CONTRAST 1: 2 -DPT. R/L/EIN R BLEND FIX 4 RO 2 R	Test no. 4	Fourth programme step
	Optotype: RO	Landolt ring gap at top right
	CONTRAST 1: 2	Presented contrast level 1:2
	-DPT. R/L/BIN: R	Monocular examination of the right eye (option), without minus lens
	BLEND FIX: (empty)	Glare and fixation aids are switched off

9.2 Filling In the Test Sheet

➔ Prepare a test sheet before every examination.

8	Name _____		Date of birth _____		1				
	First Name _____		Sex m ☐ f ☐						
	Street _____								
	City _____								
7	Examination		wears no visual aid ☐	tinted lenses ☐	binocular ☐	2			
			with spectacles ☐	clear lenses ☐	monoc. right ☐				
			with contact lenses ☐		monoc. left ☐				
6	Test No.	1	2	3	4	5	6	7	8
	Contrast	1:23	1:5	1:2,7	1:2	1:23	1:5	1:2,7	1:2
	without glare				with glare				
	at least 3 of 5 positions are recognized	☐	☐	☐	☐	☐	☐	☐	☐
5	Visual performance improves with different lense powers		-0.5 D ☐	-1.0 D ☐	-1.5 D ☐	no improvement ☐			
	Qualified for Night Driving		Yes ☐	No ☐					
4	City, date _____		Examiner _____		3				
	Mesotest II Typ 32800		for testing mesopic vision and glare sensitivity		Test form sheet No. 62855		 Optikgeräte GmbH D-35549 Wetzlar Tel. ++49/641/2005-0 Fax ++49/641/2005-55		

- | | |
|---|---|
| 1 Patient data | 5 Assessment of the night driving ability |
| 2 Examined eyes | 6 Use of concave lenses |
| 3 Information about the place, date of the examination and the person conducting the test | 7 Results of the individual tests |
| 4 Information about the device used for the examination | 8 Information about any visual aids used |
- Fig. 9-2: Test sheet

If the patient wears glasses or contact lenses for distance vision, the examination is conducted with the glasses or contact lenses.

➔ Mark this on the test sheet (including tinting, etc.), [fig. 9-2, page 25](#).

9.2.1 Starting the Examination

➔ New start: Switch on the unit.

After automatic calibration of the Mesotest II, test no. 1 is presented.

➔ After an examination: Press this button.

The programme begins with test no. 1



Instruct the patient

➔ Instruct the patient with the help of test no. 1 (contrast level 1:23), as this test is the easiest to recognise.

9.3 Test Mesopic Vision



- ➔ Press this button.
Test no. 2 is requested.

- ➔ Present five different Landolt ring positions to the patient.
You can select these randomly with the Landolt ring buttons.

The visual requirement is deemed met when the patient has recognised at least 3 of 5 positions of the Landolt ring (60% criterion).

- ➔ Mark this down on the test sheet accordingly, see [fig. 9-2, page 25](#).
- ➔ Conduct the tests, nos. 3 and 4 in the same way.

9.4 Test Mesopic Vision With Glare

When you proceed to test steps 5-8, the glare is automatically switched on in order to test the patient's mesopic vision with glare.



If the patient looks directly into the glare light when you switch over to the test for mesopic vision with glare, he may no longer be able to recognise any optotypes to begin with. This can lead to a delay in the examination procedure, as he needs some time (approx. 10-30 seconds) to readapt.

To avoid an adaptation problem, proceed as follows:



- ➔ Before you switch to test no. 5: Press this button to switch on the fixation aid.
- ➔ Tell the patient to focus on the red fixation mark. Explain to the patient that he must keep looking in that direction and must not, under any circumstances, look directly at the glare light during any of the subsequent tests.
- ➔ Press this button again to switch off the fixation aid.
Or, the fixation aids are automatically switched off for the next test or optotype, to prevent falsification of the examination results.




- ➔ Press this button.
- ➔ Conduct the remaining tests, nos. 5 through 8 in the same way as described above, [sec. 9.3, page 26](#).
- ➔ Fill out all of the other required items on the test sheet.

9.4.1 Test for Night Myopia (optional)

If your Mesotest II is equipped with minus lenses (Mesotest II konkav), you can conduct these additional examinations. With the minus lenses, you can determine whether the patient suffers from night myopia.

→ Swing in the minus lenses for the last test number that the patient just recognised.

→ Press this button to swing in a minus lens.

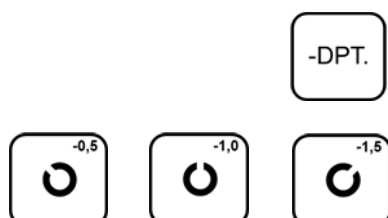

→ Press the button with the desired lens power.

If the visual acuity of the patient improves, he may suffer from night myopia. This should be corrected with suitable visual aids.

You can determine the extent of the night myopia by swinging in lenses of different powers (-0.5/-1.0/-1.5 D).

9.4.2 Monocular Examinations (optional)

If your Mesotest II is equipped with an eye occluder (Mesotest II konkav), you can conduct the following additional examinations. Monocular examinations can be conducted with the eye occluder.

This is very important, for example for patients with intraocular lenses or who have undergone excimer laser surgery, to assess the visual function of one eye.

→ Press this button to swing in the eye occluder.

→ Press the button again to advance the cover.



The information (R / L / BIN) shown in the display relates to the eye that is to be examined.

10 Conducting the Examination with a Computer/Laptop

Here, you will find information about:

- How to use the examination menu, [sec. 10.1, page 28](#)
- How to change settings, if necessary, [sec. 10.2, page 31](#)
- How to conduct an examination, [sec. 10.3, page 36](#)

10.1 Using the Examination Menu

The general use of the examination menu is explained in this chapter.

- ➔ Start the Patient Data Management, [sec. 10.3.1, page 37](#).
- ➔ If necessary, press the [Mesotest II] button to start the Mesotest II programme.
- ➔ Press the button [Examination] (1) to start an examination.


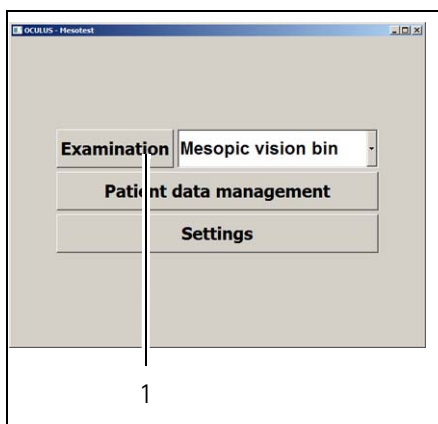


Fig. 10-1: Menu "Mesotest II"

The following components are available:

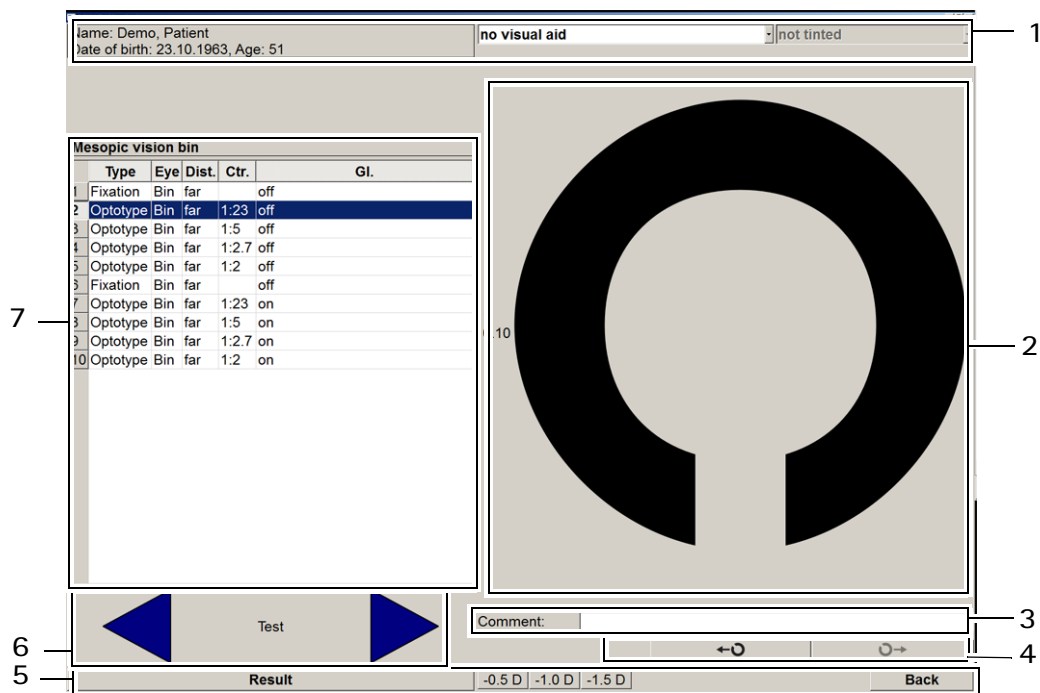

- 1 Header with information about the patient and "vision aid" selection field.
- 2 Optotype display
- 3 Input field "Comment"
- 4 Button: Change Landolt ring position
- 5 Task bar with the buttons [Back] and [Result], as well as minus lenses for myopia test (Mesotest II konkav)
- 6 Arrow keys for steps: "back" and "next"
- 7 Display the exam programme with test steps

Fig. 10-2: Examination menu

10.1.1 View the Patient's Information and Enter any Visual Aid

The name, date of birth and age of the patient are displayed in the header (item 1). If the patient wears visual aids, you must select these here before conducting the examination.

A drop-down list appears with the following options.

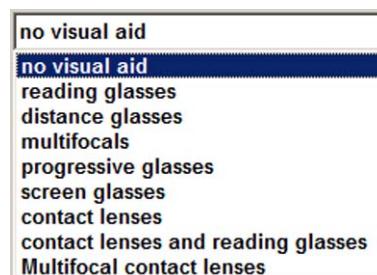


Fig. 10-3: Selection of visual aids

- ➔ Select the visual aid(s) worn by the patient.
- ➔ Specify whether the visual aid is tinted or not.

When you select a visual aid at the beginning of an examination, it is applied to all subsequent test steps.

If you select a different visual aid in a later test step, then that is applied to all subsequent test steps.

10.1.2 Marking the Results in the Test Image Display

The test image (optotype) that the patient is presently seeing is shown to you in the test image display (*fig. 10-2, page 29, item 2*). Here, you have to mark the individual test results so that the examination can be subsequently evaluated. Here, you can mark which optotypes the patient has recognised or not recognised. The following shows you how to mark the results.

- ➔ Click on the optotype if the patient has recognised it.
The recognised optotype is marked with a green tick.
If the patient does not recognise the optotype, click twice on the optotype.



Recognised optotype



Not recognised optotype

Fig. 10-4: Marking the results

10.1.3 Present the Next Optotype

1
0.10

After you have marked the results, you must present another optotype to the patient. A number beside the optotype (here: 1) indicates the first position of the optotype. You can present 5 different Landolt ring positions per contrast level to the patient, see *fig. 10-2, page 29, item 2*. The visual requirement is deemed met when the patient has recognised at least 3 of 5 positions of the Landolt ring (60% criterion).



- ➔ Press this button to present the next optotype.



- ➔ Press this button to display the previous, marked optotype.

10.1.4 Swing in Minus Lenses

You can swing in minus lenses to determine whether the patient suffers from night myopia.



You can only use the myopia test if your Mesotest II has this function (Mesotest II konkav), and the myopia test has been activated in the "Settings", see [sec. 10.1.4, page 31](#).



Fig. 10-5: The buttons for minus lenses are activated

- ➔ Swing in the minus lenses for the last test that the patient just recognised.
- ➔ Press the button with the desired lens power.

If the visual acuity of the patient improves, he may suffer from night myopia.

You can determine the extent of the night myopia by swinging in lenses of different powers (-0.5/-1.0/-1.5 D).

10.2 Change Settings

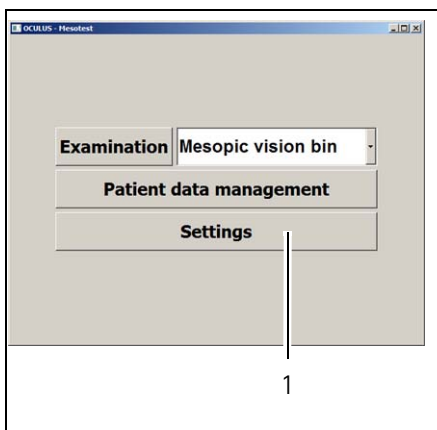


Fig. 10-6: Menu "Mesotest II"

If necessary you can make some basic adjustments in the Mesotest II programme (1). If you do not want to change any settings, proceed as described in [sec. 10.3, page 36](#).

- ➔ Select a patient from the list of patients and double-click on the selected patient to start the Mesotest II programme. The „Mesotest II“ menu is displayed.


The following settings can be changed:

- Set the date format, [page 32](#)
- Set the language, [page 32](#)
- Change the visual acuity data, [page 32](#)
- Activate the myopia test, [page 33](#)
- Enable the option to allow the patient's gender to be entered in the patient data management, [page 33](#)
- Enable full screen mode when the programme is started, [page 34](#)
- Enable to load an existing examination without using the OCULUS patient data management, [page 33](#)
- Activate Demo mode, [page 34](#)
- Edit examination programme, [page 34](#)
- Personalise the results printout, [page 34](#)

Proceed as follows:

- ➔ Press the button [Settings].

The following dialogue is displayed:

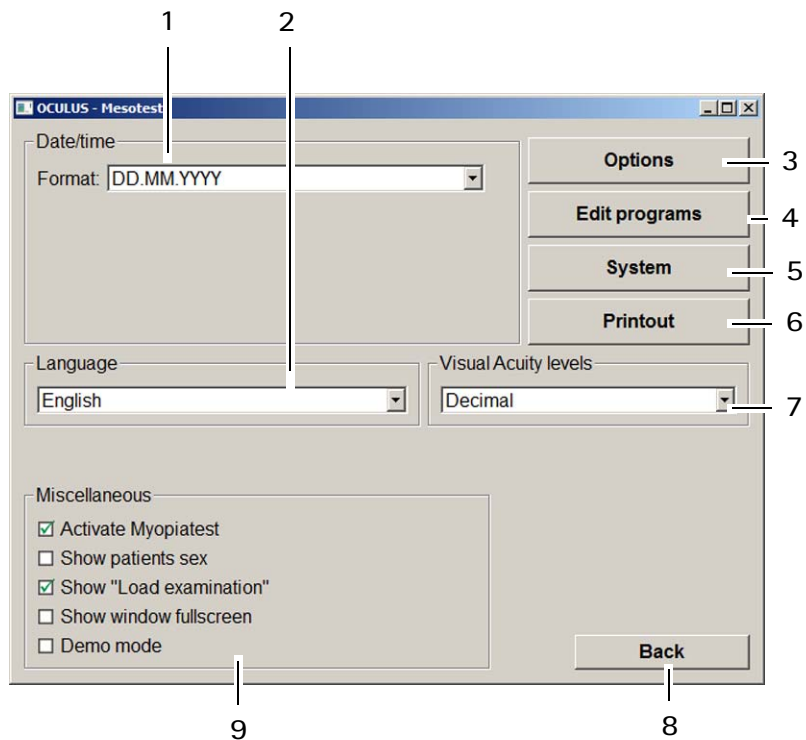

- | | |
|--------------------------|----------------------------|
| 1 "Date / Time" box | 6 [Printout] button |
| 2 "Language" box | 7 Box „Visual acuity data“ |
| 3 [Options] button | 8 [Back] button |
| 4 Button [Edit programs] | 9 "Miscellaneous" box |
| 5 Button [System] | |

Fig. 10-7: "Settings" screen

Change date format

"Date / Time" box (fig. 10-7, page 32, item 1):

- ➔ From the drop-down list, select the desired date format.

Language settings

"Language" box (fig. 10-7, page 32, item 2):

- ➔ Select the desired language from the drop-down list.
The selected language will be active after the next start of the Mesotest II programme.

Change visual acuity data

"Visual acuity levels" box (fig. 10-7, page 32, item 7):

- ➔ Select the desired setting for the visual acuity data from the drop-down list.

Activate myopia test

"Miscellaneous" box (fig. 10-7, page 32, item 9):

This test is useful for determining whether the patient suffers from night myopia.



You can only use this test if your Mesotest II has this function (Mesotest II konkav).



Fig. 10-8: The buttons for minus lenses are activated

- ➔ Activate the checkbox [Activate Myopiatest].
Three additional buttons,
-0,5 D / -1,0 D / -1,5 D, now appear on the user interface.

Show patients sex

"Miscellaneous" box (fig. 10-7, page 32, item 9):

- ➔ Activate the [Show patients sex] checkbox.
The patient's gender is displayed in the examination menu.
Prerequisite: The feature „Patient's sex" was entered and activated in the Patient Data Management (sec. 12.5, page 57).

Display: Show „Load examination"

If you do not work with the OCULUS patient data management, you can use this function to load an existing examination.

"Miscellaneous" box (fig. 10-7, page 32, item 9):

- ➔ Enable the [Show "Load examination"] checkbox.
The according button (1) will added in the "Mesotest II" menu.

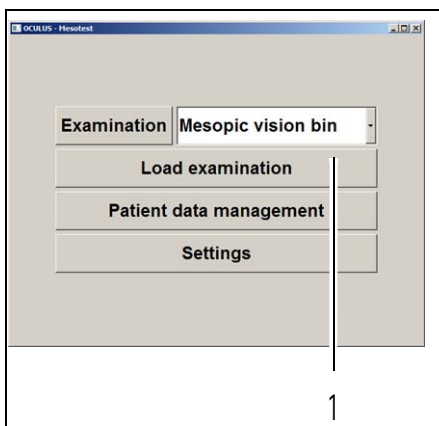


Fig. 10-9: [Load examination] button

Show window fullscreen

"Miscellaneous" box ([fig. 10-7, page 32, item 9](#)):

- ➔ Activate the [Show window fullscreen] checkbox.

Activate "Demo" mode

This function is for demonstration purposes only, if a Mesotest II is not connected to your laptop/netbook/PC.

"Miscellaneous" box ([fig. 10-7, page 32, item 9](#)):

- ➔ Activate the [Demo mode] checkbox.

Edit examination programme

- ➔ Press the button [Edit programs] ([fig. 10-7, page 32, item 5](#)).

You can

- change the order in which the examination programmes are displayed
- copy and edit an existing examination programme
- create and edit a new examination programme


You will find more detailed information in [sec. 11, page 47](#).

Personalise the results printout

- ➔ Press the button [Result] ([fig. 10-7, page 32, item 5](#)). On this screen, you can enter data about yourself or your office address, and add your own logo to the results printout.

Screen options

This screen is for information purposes.

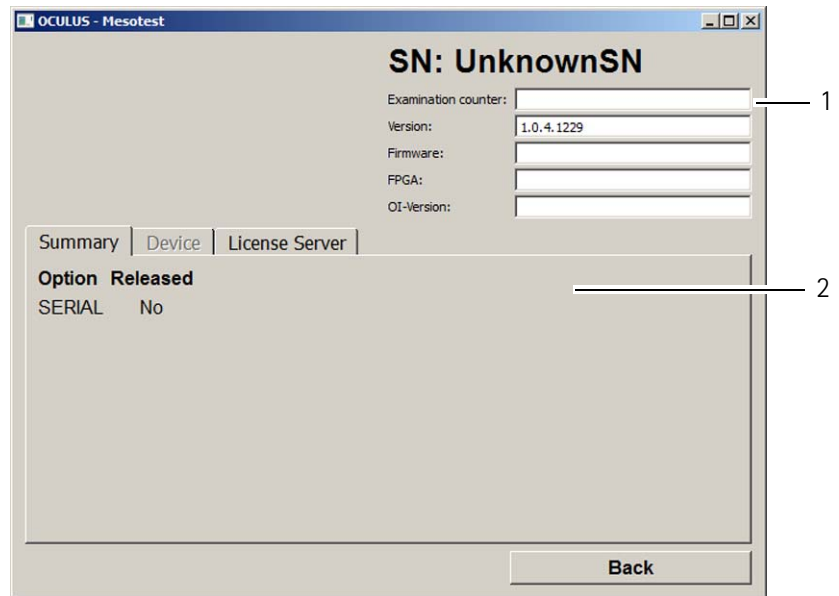

1 Device information

2 Connection status to the floating license key

Fig. 10-10: "Options" screen

System screen

You can activate the interface to external office management programmes here.

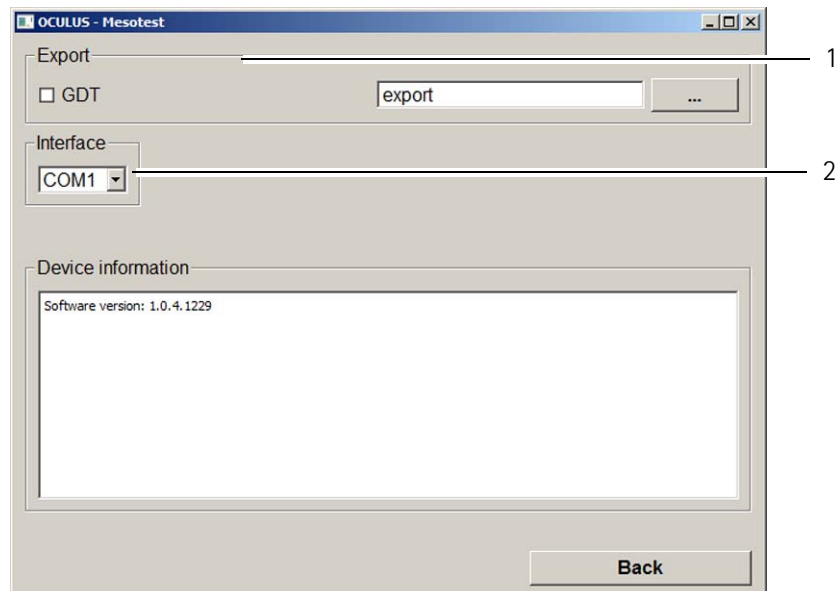


Fig. 10-11: „System“ screen

- ➔ Activate the checkbox [GDT] in the „Export“ box (1).
- ➔ Press the button [...] in the „Export“ box (1).

- ➔ In the dialogue box, select the location to which the data are to be saved.

This interface (2) is used to establish a connection between the Mesotest II and the computer.

- ➔ Select the appropriate interface port.
Refer to the Operating Manual for your computer to find your computer's COM port numbers.
A notebook generally has only one port: COM1. A desktop computer generally has two: COM1 and COM2.
- ➔ Save your selection.

Switch to the Mesotest II programme

- ➔ Press the [Back] button (*fig. 10-7, page 32, item 8*).
The settings will take effect immediately, except for the language change.

10.3 Conduct the Examination

To conduct an examination:

- Start the Patient Data Management, *sec. 10.3.1, page 37*
- Select and start an examination, *sec. 10.4, page 39*
- Test mesopic vision, *sec. , page 41*
- Test mesopic vision with glare, if applicable, *sec. , page 41*

You can display and print out the results, *sec. 10.5, page 45*.

- ➔ Prepare the patient for the examination, *sec. 8, page 20*.

10.3.1 Starting Patient Data Management


You can enter patient data in the Patient Data Management and then use it.

After you have switched on the computer, it first loads the operating system.

If the Windows desktop is displayed, you have to start the Patient Data Management programme from there.

➔ If necessary, click on the Mesotest II icon: .

The user interface for the Patient Data Management is displayed

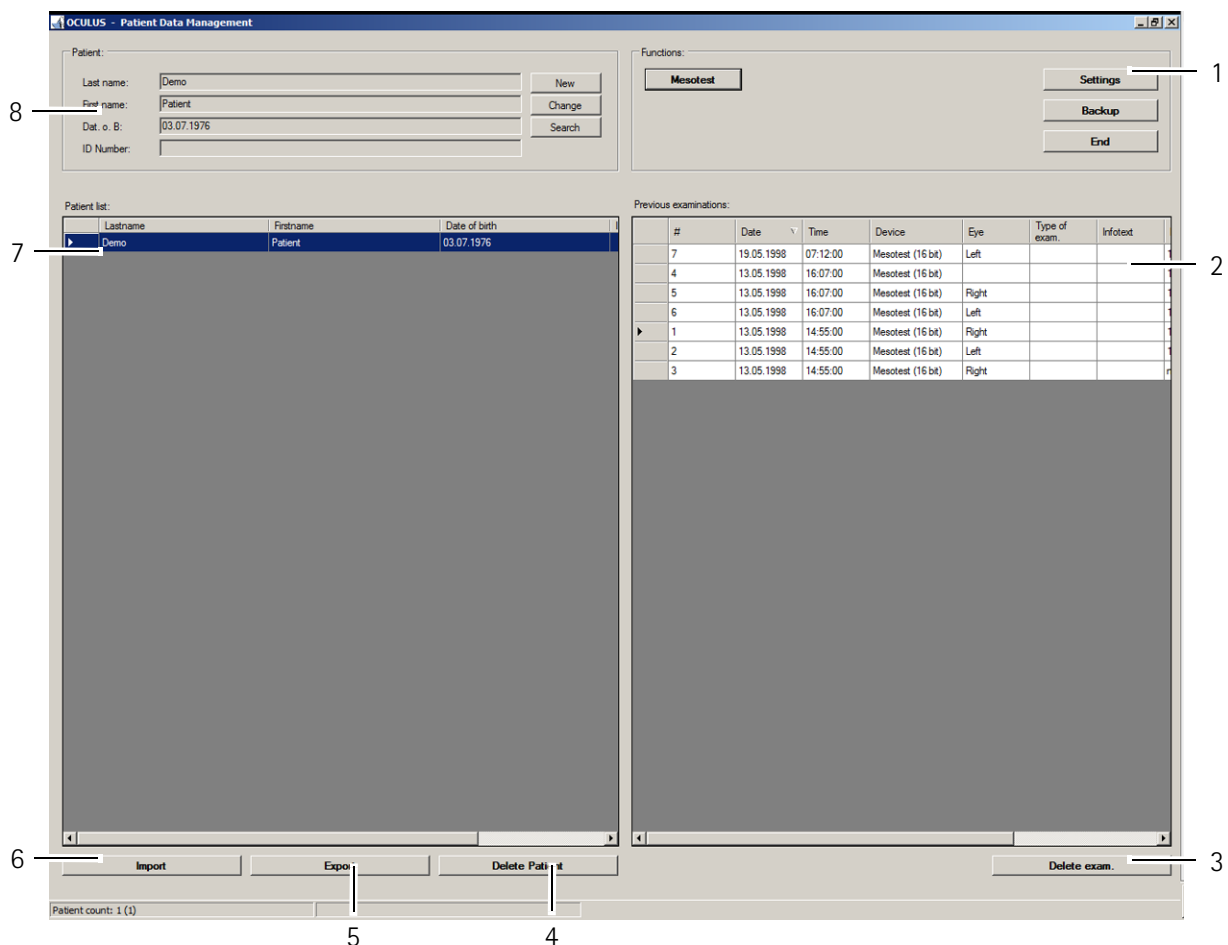

- 1 "Functions" box
- 2 Previous examinations
- 3 [Delete exam] button
- 4 [Del. Patient] button

- 5 [Export] button
- 6 [Import] button
- 7 Patient list
- 8 "Patient" box


Fig. 10-12: User Interface for Patient Data Management



To start the Mesotest II programme, first enter a new patient (8), or select an existing patient from the list of patients (7). For more information on Patient Data Management, refer to the [sec. 11, page 47](#).

10.3.2 Entering a New Patient

- ➔ Press the [New] button to enter a new patient in the Patient Data Management system.
- ➔ Enter the patient's last name, first name and date of birth in the patient window.

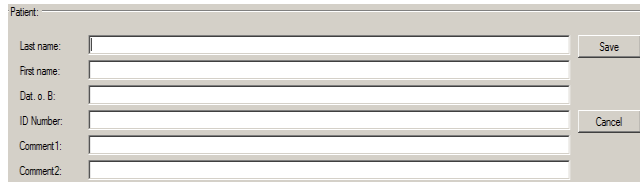

The form is titled 'Patient:' and contains several input fields and two buttons. The fields are: 'Last name:', 'First name:', 'Dat. o. B:', 'ID Number:', 'Comment1:', and 'Comment2:'. The 'Save' button is located to the right of the 'Last name' field, and the 'Cancel' button is located to the right of the 'ID Number' field.


Fig. 10-13: Entering patients

Optionally you can enter an ID number for the patient.

- ➔ To save the data you entered, click [Save].
- The patient you have just entered now appears in the patient list.

10.3.3 Selecting an Existing Patient

The patient data list on the left-hand side of the screen displays all previously examined patients in alphabetical order.

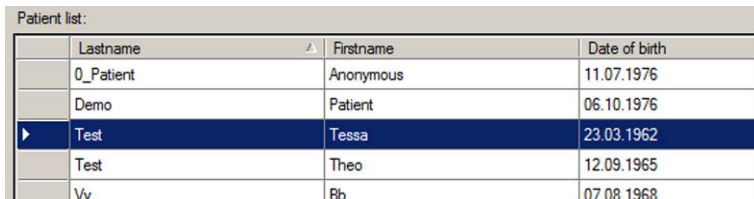

Patient list:			
	Lastname	Firstname	Date of birth
	0_Patient	Anonymous	11.07.1976
	Demo	Patient	06.10.1976
▶	Test	Tessa	23.03.1962
	Test	Theo	12.09.1965
	Vv	Rh	07.08.1968

Fig. 10-14: Patient list

- ➔ Choose [Search] to quickly find the patient you are looking for in the list.
 - ➔ Enter the patient's name or the first letter of the name in the "Last name" field.
- Alternatively, you can search for the patient using an ID number, first name or date of birth, assuming that one was assigned when the patient was first recorded.
- ➔ To transfer the patient's name to the patient window, click the entry that you need in the list. This also brings up a list of any previous examinations for that patient in the examination window (bottom right side).

Extended patient search: [Extended] checkbox

- ➔ Click on the [Extended] checkbox.

The screen displays additional search parameters which reference previous examinations.

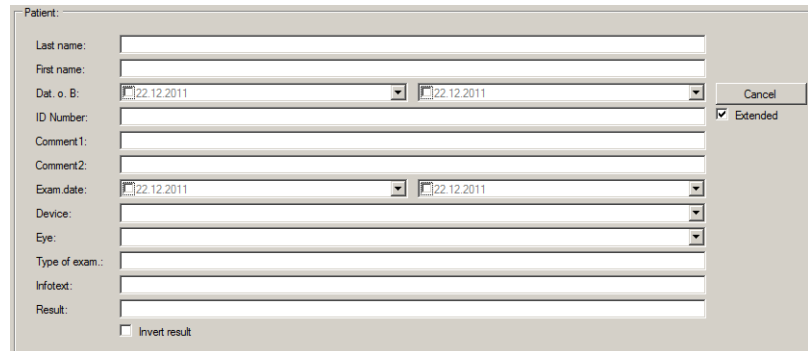


Fig. 10-15: Advanced search

- ➔ Enter a search parameter.
The respective patient is highlighted in the list.
- ➔ Press the button [End search].
- ➔ Proceed as for entering a patient name.

10.4 Select and Start Exam

- ➔ If necessary, press the [Mesotest II] button to start the Mesotest II programme.

The following dialogue appears:

- ➔ Select a program (2):
You can choose between "Mesopic vision bin" and „Mesopic vision mono“.

If you have created your own programs, a drop-down list showing all programs is displayed here, [sec. 11, page 47](#).

- ➔ Press the button [Examination] (1) to start the examination.

The selected program appears:

- Mesopic vision bin, [sec. 10.4.1, page 40](#)
- Mesopic vision mono (optional with Mesotest II konkav), [sec. 10.4.2, page 42](#)


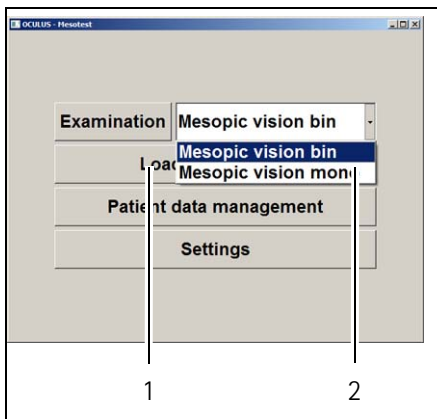


Fig. 10-16: Menu "Mesotest II"

10.4.1 Conduct the Programm "Mesopic vision bin"

With the examination program „Mesopic vision bin“, you can examine mesopic vision binocularly, at different contrast levels, with and without glare. It consists of 10 test steps.

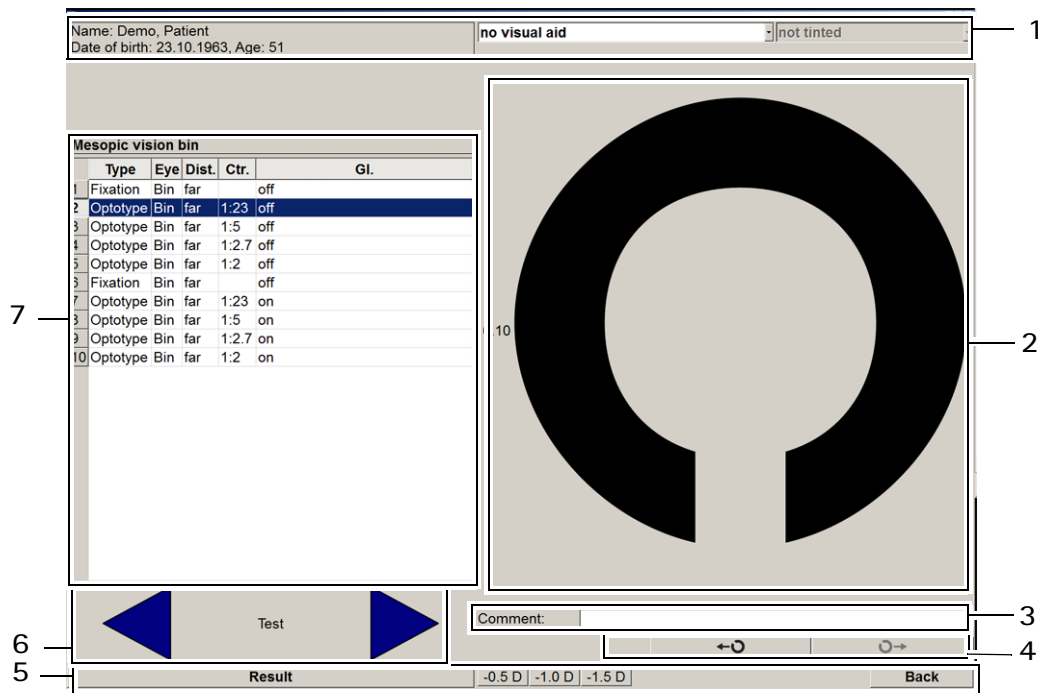

- 1 Header with information about the patient and "vision aid" selection field.
- 2 Optotype display
- 3 Input field "Comment"
- 4 Button: Change Landolt ring position
- 5 Task bar with the buttons [Back] and [Result], as well as minus lenses for myopia test (Mesotest II konkav)
- 6 Arrow keys for steps: "back" and "next"
- 7 Display the exam programme with test steps

Fig. 10-17: Examination menu

Test 1 and Test 6: Facilitate the fixation, as here the fixation marks are displayed

Test 2-5: Are used to test mesopic vision

Test 7-10: Are used to test mesopic vision with glare.

Enter a visual aid

If the patient wears glasses or contact lenses, the examination is conducted with the glasses or contact lenses.

- ➔ Enter whether the patient wears visual aids (including any tinting, etc.), [sec. 10.1.1, page 29](#).

Instruct the patient

- ➔ Select test step 1

- Tell the patient to focus between the red fixation marks.
- Instruct the patient with the help of displayed Landolt ring.

Test mesopic vision

- Select test step 2.
- Present the optotype to the patient in 5 different positions, [sec. 10.1.3, page 30](#).
- Mark the results, [sec. 10.1.2, page 30](#).

The visual requirement is deemed met when the patient has recognised at least 3 of 5 positions of the Landolt ring (60% criterion).

- Conduct test steps 3 through 5 in the same way.


Test mesopic vision with glare

- Select test step 6.
- Tell the patient to focus on the red fixation mark. Explain to the patient that he must keep looking in that direction and must not, under any circumstances, look directly at the glare light during any of the subsequent tests.
- Conduct test steps 7 through 10 in the same way as described above, [sec. , page 41](#).

The visual requirement is deemed met when the patient has recognised at least 3 of 5 positions of the Landolt ring (60% criterion).

10.4.2 Conduct the Programm "Mesopic vision mono"

With the examination program „Mesopic vision mono“ (optional with Mesotest II konkav), you can examine mesopic vision monocularly, for the right and for the left eye, at different contrast levels, with and without glare. It consists of 20 test steps.

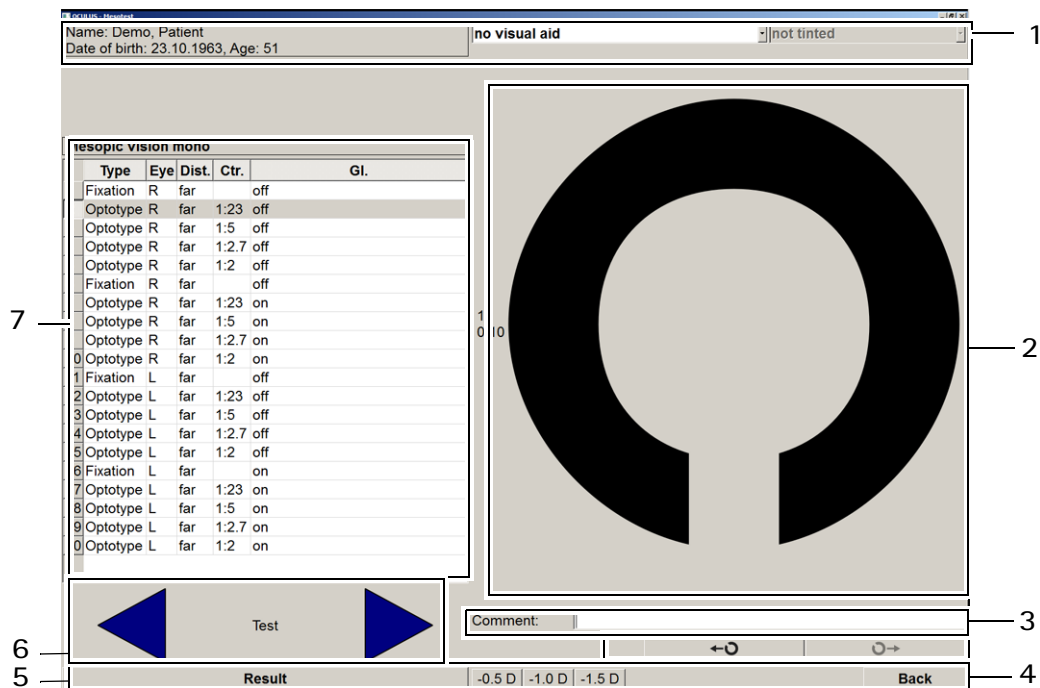

- 1 Header with information about the patient and "vision aid" selection field.
- 2 Optotype display
- 3 Input field "Comment"
- 4 Button: Change Landolt ring position
- 5 Task bar with the buttons [Back] and [Result], as well as minus lenses for myopia test (Mesotest II konkav)
- 6 Arrow keys for steps: "back" and "next"
- 7 Display the exam programme with test steps

Fig. 10-18: Examination menu

Test 1, 6, 11 and 16: Facilitate the fixation, as here the fixation marks are displayed

Test 2-5 and 12-15: Are used to test mesopic vision

Test 7-10 and 17-20: Are used to test mesopic vision with glare.

Enter a visual aid

If the patient wears glasses or contact lenses, the examination is conducted with the glasses or contact lenses.

➔ Enter whether the patient wears visual aids (including any tinting, etc.), [sec. 10.1.1, page 29](#).

Instruct the patient

- ➔ Select test step 1
- ➔ Tell the patient to focus between the red fixation marks.
- ➔ Instruct the patient with the help of displayed Landolt ring.

Test mesopic vision for the right eye

- ➔ Select test step 2.
- ➔ Present the optotype to the patient in 5 different positions, [sec. 10.1.3, page 30](#).
- ➔ Mark the results, [sec. 10.1.2, page 30](#).

The visual requirement is deemed met when the patient has recognised at least 3 of 5 positions of the Landolt ring (60% criterion).

- ➔ Conduct test steps 3 through 5 in the same way.

Test mesopic vision with glare for the right eye

- ➔ Select test step 6.
- ➔ Tell the patient to focus on the red fixation mark. Explain to the patient that he must keep looking in that direction and must not, under any circumstances, look directly at the glare light during any of the subsequent tests.
- ➔ Conduct test steps 7 through 10 in the same way as described above, [sec. , page 41](#).

The visual requirement is deemed met when the patient has recognised at least 3 of 5 positions of the Landolt ring (60% criterion).

Test mesopic vision for the left eye

- ➔ Select test step 11.
- ➔ Tell the patient to focus between the red fixation marks.
- ➔ Present the optotype to the patient in 5 different positions (step 12), [sec. 10.1.3, page 30](#).
- ➔ Mark the results, [sec. 10.1.2, page 30](#).

The visual requirement is deemed met when the patient has recognised at least 3 of 5 positions of the Landolt ring (60% criterion).

- ➔ Conduct test steps 13 through 15 in the same way.

Test mesopic vision with glare for the left eye

- ➔ Select test step 16.
- ➔ Tell the patient to focus between the red fixation marks. Explain to the patient that he must keep looking in that direction and must not, under any circumstances, look directly at the glare light during any of the subsequent tests.
- ➔ Conduct test steps 17 through 20 in the same way as described above, [sec. , page 41](#).

The visual requirement is deemed met when the patient has recognised at least 3 of 5 positions of the Landolt ring (60% criterion).

10.4.3 Test for Night Myopia (optional)

If your Mesotest II is equipped with minus lenses (Mesotest II, konkav), you can conduct these additional examinations. With the minus lenses, you can determine whether the patient suffers from night myopia.

The myopia test must be activated, "[Activate myopia test" auf Seite 33](#).

- ➔ Swing in the minus lenses for the last test that the patient just recognised.



Fig. 10-19: The buttons for minus lenses are activated

- ➔ Press the button for the desired lens power.

If the visual acuity of the patient improves, he may suffer from night myopia.

You can determine the extent of the night myopia by swinging in lenses of different powers (-0.5/-1.0/-1.5 D).

10.5 Display and Print Results

Result

You can view the test results.

➔ To do so, press this button.



Then you can print the test results.

➔ To do so, press this button.

OCULUS Mesotest

OCULUS Optikgeräte GmbH
Münchholzhäuserstr. 29
35582 Wetzlar



Patient: Mustermann, Franz (M)		Date of exam.: 01.08.2014	
Date of birth: 08.08.1962		Program: Dämm. Blend. Bin	
Age: 51		Comment:	

1

Contrast level	Glare	Seen	Distance	Visual aid	Concav lenses	Eye
1:2.7	without	3 of 5 (60%)	far	progressive glasses (not tinted)		Both
1:2.7	with	3 of 5 (60%)	far	progressive glasses (not tinted)		Both

2

Contrast level	Glare	Seen	Distance	Visual aid	Concav lenses	Eye
	without	with addition of concav lenses no higher level of contrast achieved				Both
	with	with addition of concav lenses no higher level of contrast achieved				Both

3

Remark:
If an additional concave lens improves the contrast vision, an indication for night myopia is given.

4

1

1 Header with information about patients

3 Test results „Mesopic vision and glare testing“ when testing was done using a Mesotest II konkav

2 Test results „Mesopic vision and glare testing“

4 Remark

Fig. 10-20: Results printout

➔ To return to the examination: Press the „Close“ X

10.6 Finishing the Exam

- ➔ Press the button [Back] to end the examination.
 - ➔ In the dialogue, select whether the data are to be saved or not.
- You are then returned to the Mesotest II main menu. You can now
- Start a new examination ([sec. 10.4, page 39](#)).
 - Change patients. To do this, press the button [Select patient].
 - ➔ Clean and disinfect viewer enclosure after each exam, [sec. 13.2, page 69](#).

10.7 Load Existing Examination

You can load an existing examination, e.g. in order to print it out.

If you are working with the OCULUS Patient Data Management, you can load an existing examination from there.

- ➔ Double click on the desired examination in the list of examinations, [fig. 10-12, page 37, item 2](#).

The examination is displayed in the examination menu

If you are not working with the OCULUS Patient Data Management, you can nevertheless load an existing examination. To do so, you must have activated the appropriate button in the settings, ([fig. 10-7, page 32, item 9](#)).

- ➔ Press the button (1).

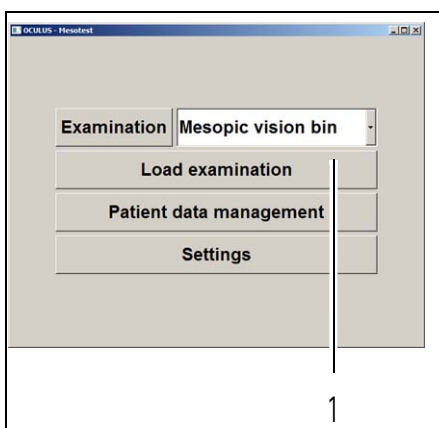


Fig. 10-21: [Load examination] button

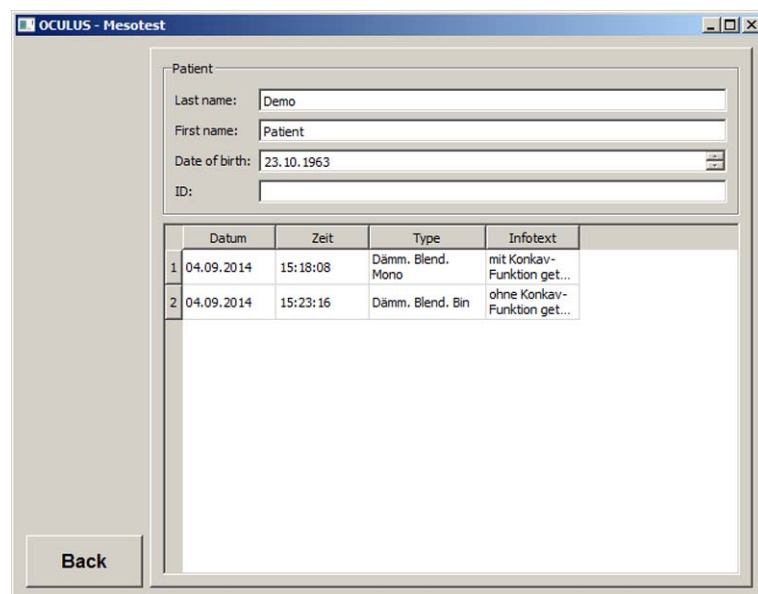


Fig. 10-22: Loading an examination without using the OCULUS Patient Data Management

- ➔ Move the cursor to the desired examination.
- ➔ Double-click on that examination.
- ➔ The examination is displayed in the examination menu.
- ➔ If you do not want to load any of the examinations, press the [Back] button. This returns you to the Mesotest II menu.

11 Editing the Examination Programmes

With the Mesotest II, you can

- Compile a custom selection list of examination programmes, [sec. 11.1, page 47](#)
- Create a new, custom examination programme, [sec. 11.2, page 48](#)
- Edit an examination programme that you created, [sec. 11.3, page 51](#)

11.1 Compilation of a Custom Selection List of Examinations

Mesotest II > Settings > Edit programmes

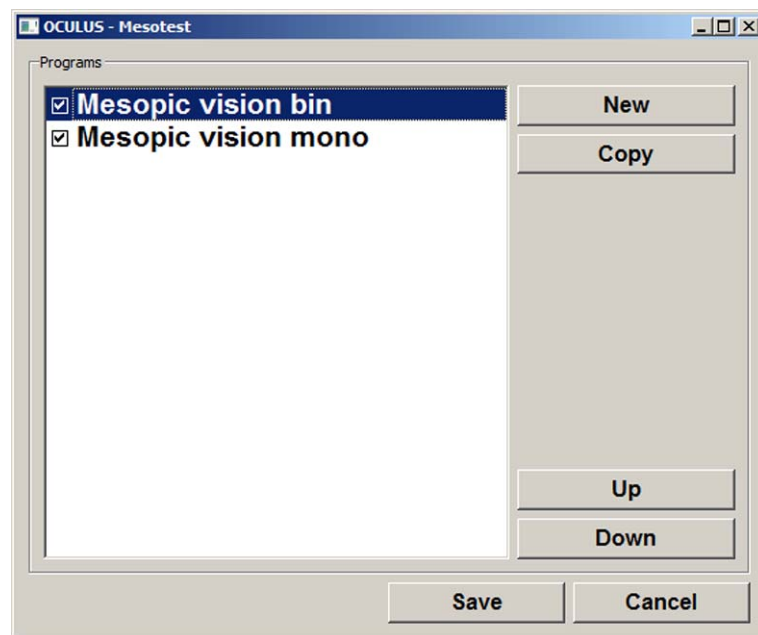


Fig. 11-1: Activation of an examination programme

This function acts on the display in the selection list for the examination programmes.

You can activate and deactivate the wanted examination programmes and can change the order in which they are displayed here.

Activate/deactivate examination programme

- ➔ Activate/deactivate the checkbox for the wanted examination programme.
- ➔ Press the [Save] button.
The selection list is adapted accordingly.

Change order

- ➔ Select an examination
- ➔ Press the [Up] or [Down] button as many times as necessary until the examination has moved to the wanted position.
- ➔ Press the [Save] button.
The selection list is adapted accordingly.

11.2 Create A New Examination Programme

Mesotest II > Settings > Edit programmes

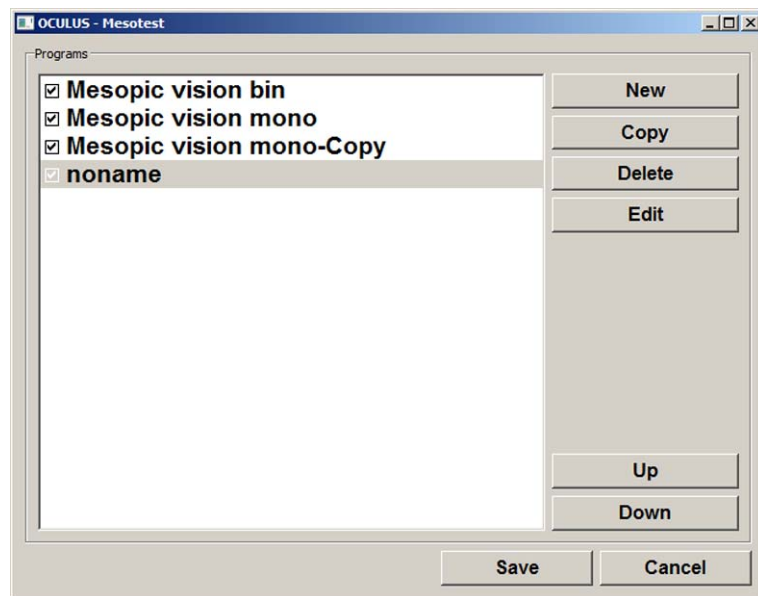


Fig. 11-2: List of activated/deactivated examination programmes.

- ➔ Press the [New] button.
A new programme appears at the end of the list under the name „noname“. The [Delete] and [Edit] buttons appear.

➔ Press the [Edit] button. The following dialogue is displayed:

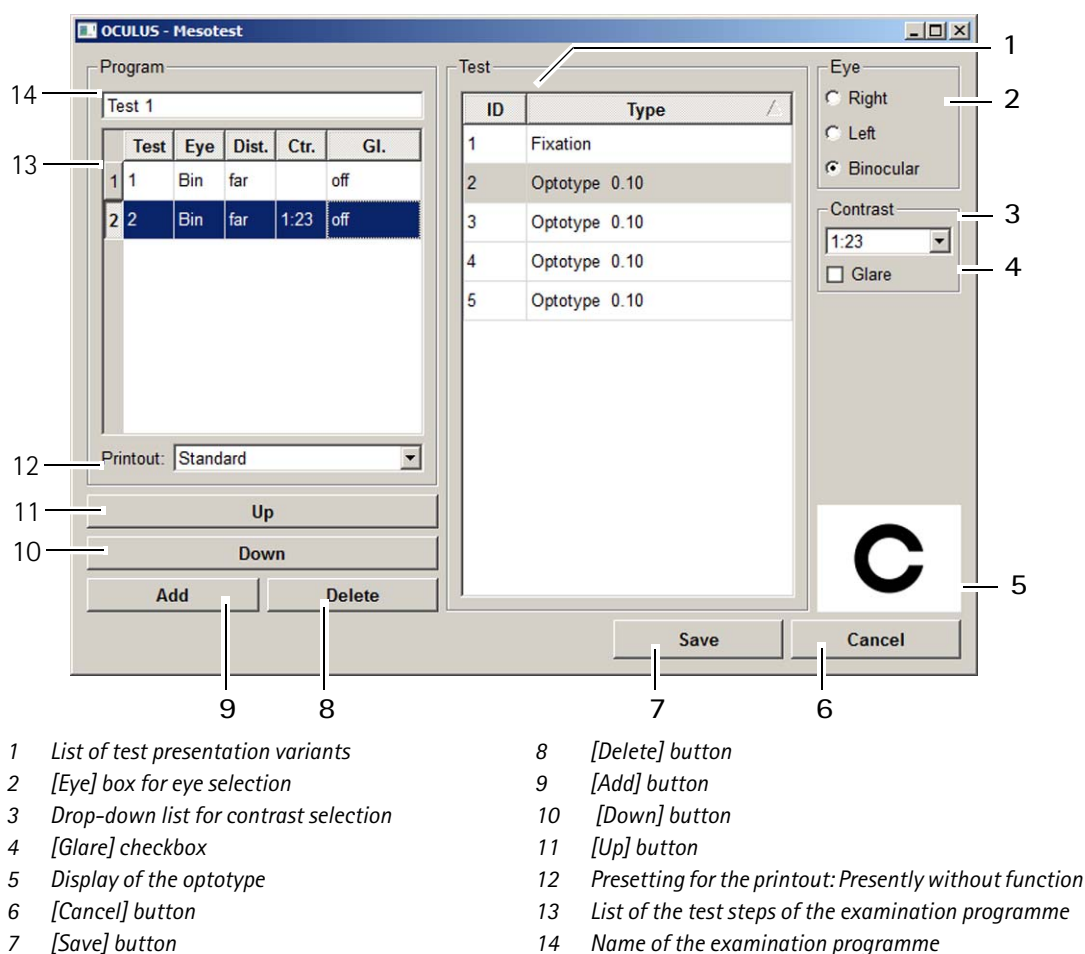


Fig. 11-3: Editing a new examination programme

Assign a name

You can assign a new name to the examination programme in the field for the "name of the examination programme" (14).

➔ To do this, move the cursor into that field and enter the new name.

Change test step settings

- ➔ In the list of test steps (13), highlight the test step for which you want to change the settings.
- ➔ Select a test presentation variant from the list (1).

Test	
ID	Type
1	Fixation
2	Optotype 0.10
3	Optotype 0.10
4	Optotype 0.10
5	Optotype 0.10

Program					
Test 1					
	Test	Eye	Dist.	Ctr.	Gl.
1	1	Bin	far		off
2	2	Bin	far	1:23	off

This is taken over into the table for the highlighted test step. The test number is shown in the first column („Test“).

Select the eye to be examined (2)

Eye
☐ Left
☐ Right
☒ Binocular

- ➔ Activate the checkbox for the eye that is to be examined (2).
- All test steps can be conducted binocularly, or at option monocularly (only possible with device type Mesotest II konkav).
- The selection is displayed in the list of test steps for the examination programme (13).

Select the contrast (3)

1:23
1:23
1:5
1:2.7
1:2

Contrast
1:23
☒ Glare

You can select any of the contrast levels shown here.

To test mesopic vision with glare:

- ➔ Activate the checkbox [Glare].

Add a test step

Test	
ID	Type
1	Fixation
2	Optotype 0.10
3	Optotype 0.10
4	Optotype 0.10
5	Optotype 0.10

To add another test step to the examination programme:

- ➔ Press the [Add] button.

You can now edit the added test step in the described manner.

Save examination programme

Program					
Test 1					
	Test	Eye	Dist.	Ctr.	Gl.
1	1	Bin	far		off
2	2	Bin	far	1:23	off

- ➔ Enter a name for the created examination programme.

- ➔ Press the button [Save] to save your examination programme.

The newly created examination programme is displayed in the list of activated and deactivated examination programmes. It is automatically set to active, [sec. 11.1, page 47](#).

11.3 Edit Examination Programme

Mesotest II > Settings > Edit programmes

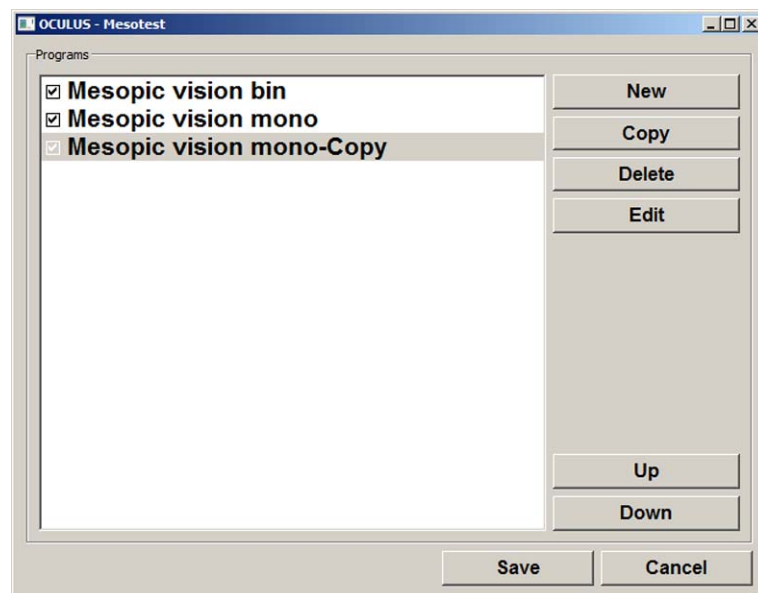


Fig. 11-4: List of activated/deactivated examination programmes.

You can only edit or delete those examination programmes that were created as user-defined ones.

- ➔ Select an examination programme.

When you select an examination programme that has already been edited or a newly created programme, the buttons [Edit] and [Delete] are displayed.

- ➔ Press the [Edit] button.
- ➔ Proceed as described in [sec. 11.2, page 48](#).

If you select a preset (default) examination programme, the buttons [Edit] and [Delete] are not displayed.

Proceed as follows.

- ➔ Press the [Copy] button.
The examination programme is copied and is displayed with the suffix "Copy" at the end of the list of programmes.
The buttons [Edit] and [Delete] appear.
- ➔ Press the [Edit] button.
- ➔ Proceed as described in [sec. 11.2, page 48](#).

12 Patient Data Management

You can do the following to the patient data

- rename, [sec. 12.1, page 52](#)
- export, [sec. 12.2, page 53](#)
- import, [sec. 12.3, page 54](#)
- save, [sec. 12.4, page 55](#)

You can also change the Patient Data Management settings, [sec. 12.5, page 57](#).

12.1 Rename Patient Data


Patient data can be changed retroactively after it has been added.

- ➔ Press the [Change] button.
The input boxes for patient data are now enabled, and the cursor jumps to the "Last name" field.
- ➔ Change the entries in the individual boxes.
- ➔ Press the [Save] button.

12.2 Exporting Patient Data

For example, patient and examination data can be exported for forwarding to another clinic.

- ➔ Select the patient and also one of the examinations in the respective list.
- ➔ Press the [Export] button below the patient list. The following dialog appears:

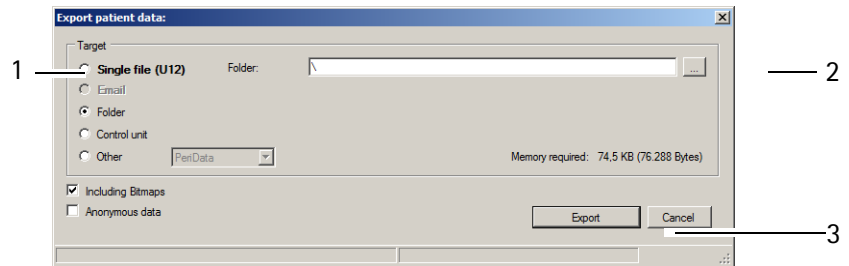

- 1 Select destination
- 2 [...] button
- 3 [Cancel] and [Export] buttons

Fig. 12-1: "Export patient data" dialog



The options for import and export of data are set as defaults in the "settings" field, see also [sec. 12.5, page 57](#).

Depending on the settings you may not have to perform all of the following steps (for example selection of the directory).

- ➔ Select the "Target" (1) where you would like to export the data.



Recommendation: Export the patient data using the "Individual file (U12)" option.

- ➔ Press the [...] button. (2)
- ➔ In the dialog that appears, select the folder or the file to which the patient data should be exported.
- ➔ Confirm your selection with [OK] or [Save].
- ➔ To export the data, press the [Export] button (3).

12.3 Importing Patient Data

If you keep patient data on a USB stick, you can import this data. Generally speaking, the various software versions are compatible with each other, so that you can still import data, even if they have been exported from a newer or an older version of the Patient Data Management system.




Note

Loss of data due to computer virus

Computer viruses can cause loss of data.

➔ Run a virus check before importing data from the USB stick.

➔ Press the button [Import]. The following dialog appears:

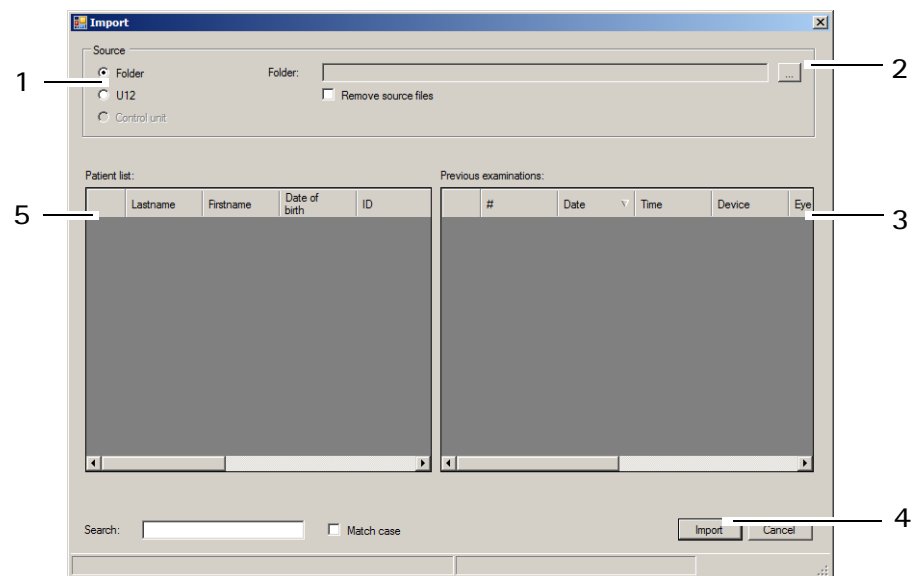

1 Select the source of the data

4 [Import] button

2 [...] button

5 Patient list

3 Previous examinations

Fig. 12-2: "Import" dialog



The options for import and export of data are set as defaults in the "settings" field, see also [sec. 12.5.3, page 62](#).

Depending on the settings you may not have to perform all the following steps (for example selection of the directory).

➔ Select the option (1) which contains the source data ("Folder" or "Single file (U12)").



Recommendation: Import the patient data using the "Single file (U12)" option.

- ➔ Press the [...] button (2).
- ➔ In the dialog box, select the directory or the file where the patient data is located.
- ➔ Confirm your selection with [OK] or [Save].
The patients that are located and the associated examinations are displayed in the lower part of the dialog.
- ➔ To import the data, press the [Import] button (4).
The data will then be available in the patient data management system.

12.4 Data Backup

You should carry out a backup of patient and examination data at regular intervals. In case of a loss of data, you can reconstruct the data from a previously created backup with the help of this function. Since data backup takes several minutes depending on the scope of the database and the data to be backed up, a backup should be carried out when the computer and the device will not be needed.



Note

Loss of data due to computer viruses

Computer viruses can cause loss of data.


- ➔ Run a virus check before importing data from the USB flash drive.



The general rules for security backups apply to backup copies created with the help of the Patient Data Management user interface. Storage of backup files should always be done on a separate system (for example on a USB flash drive with adequate capacity).

12.4.1 Backup data

- ➔ Press the [Backup] button on the upper right part in the Patient Data Management user interface. The following dialog appears:

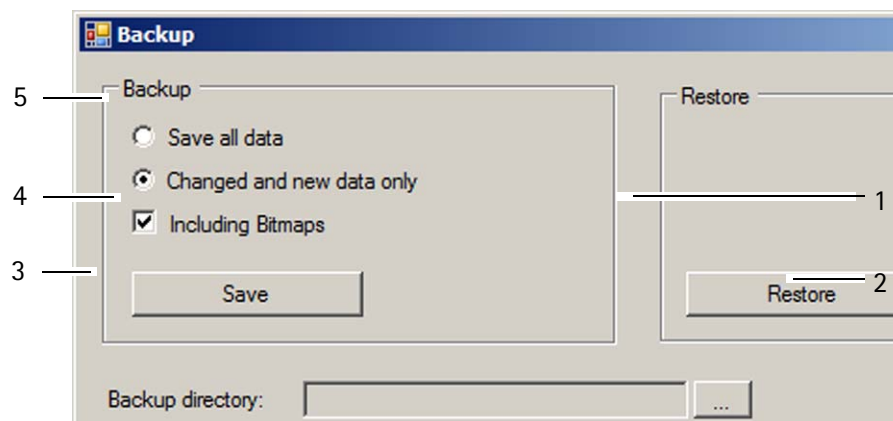

- | | |
|----------------------------------|-------------------------|
| 1 [Restore] button | 4 [Save] button |
| 2 Display free storage space. | 5 Backup data selection |
| 3 Backup folder and button [...] | |
- Fig. 12-3: "Backup" dialog

- ➔ Select whether all of the data or only changed data should be backed up.



The Patient Data Management function internally tags all saved data records.

If you selecting the option "Changed and new data only", only the data records that were not saved during a previous backup will be backed up.

- ➔ Press the [...] button to the right of the "Backup directory" box (3).
- ➔ In the dialog that appears, select the folder to which the data should be backed up.
- ➔ Confirm your selection with [OK].
- ➔ To back up the data, press the [Save] button (4). The previously selected data will then be backed up to the corresponding folder.

12.4.2 Reconstruct data

If a loss of data occurs, the data from a previous backup can be re-imported into the Patient Data Management user interface.

- ➔ Press the [...] button to the right of the "Backup directory" box (3).
- ➔ In the dialog that appears, select the folder which contains the backup data.
- ➔ Confirm your selection with [OK].
- ➔ To import the data, press the [Restore] button (1). All data in the appropriate directory are copied to the Patient Data Management software.

12.4.3 Automatic backup

In addition to the manually performed backup, it is also possible to automatically run a backup when exiting the Patient Data Management program. The settings required for this can be made in the "Settings" area, see [sec. 12.5.3, page 62](#).

12.5 Change settings

Basic specifications for working with the Patient Data Management user interface can be made in the "Settings" area.

➔ Press the [Settings] button in the upper right portion of the Patient Data Management user interface.

The "Settings" menu page will be displayed. The "Main" tab will appear in the foreground.


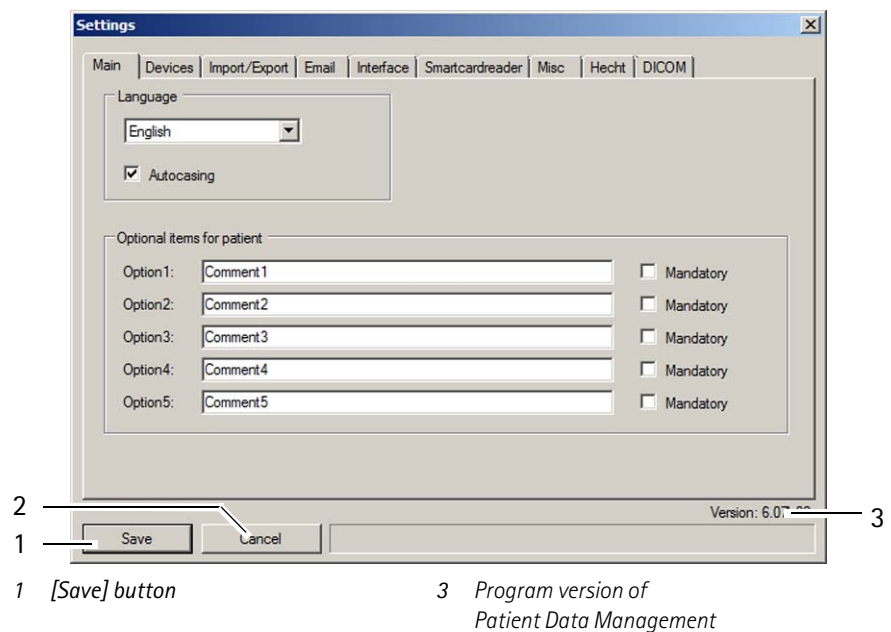


Fig. 12-4: "Settings" screen

The following information and buttons are available to you on all tabs:

- The version of the Patient Data Management program appears on the lower right (3).
- There are two buttons below and to the left for saving (1) or discarding (2) the changes that have been made. All changes are always saved or discarded, and then the screen is closed.

12.5.1 "Main" tab

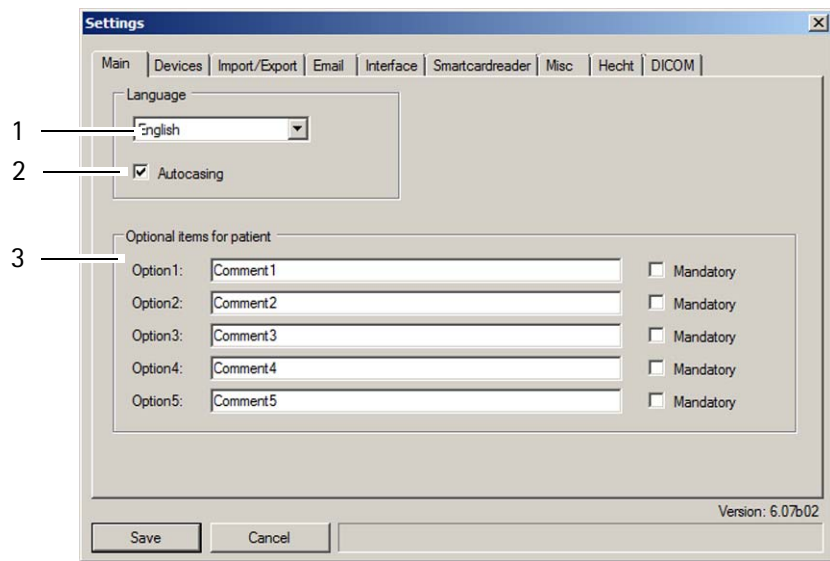

- 1 "Language selection" drop-down list 3 Optional patient entries
 2 Checkbox [automatic upper case /lower case]

Fig. 12-5: "Settings" screen, "Main" tab

"Language" box

- ➔ Select the language to be used to display the Patient Data Management user interface from the "Language" drop-down list (1).
- ➔ Enable the [Autocasing] checkbox (2) as required. If the checkbox is active, the first letter of a patient's first and last name are **always** converted to capital letters.

"Optional items for patient" box


Besides the five standard attributes of first name, last name, date of birth, gender, and ID, up to five additional attributes can be freely defined.

- ➔ Enter the identifiers for the attributes in the fields Option 1 to 5, for example, "Comments".

To be able to carry out entries for the newly defined attributes, proceed as follows:

- ➔ Click with the right mouse button on the patient list and open the associated context menu.

➔ Select the desired attribute, for example "Comment2" [2].

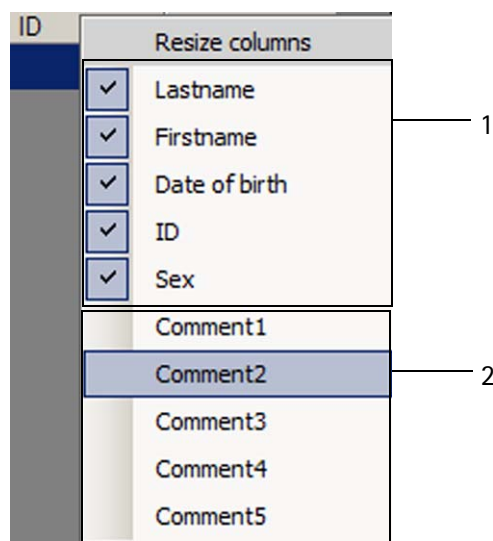


1 Previously enabled attributes

2 New attribute selected

Fig. 12-6: New attribute enabled

The context menu will close, and the attribute "Comment2" will be displayed in the upper section of the entry fields for patient data (1) as well as in the patient list (2).

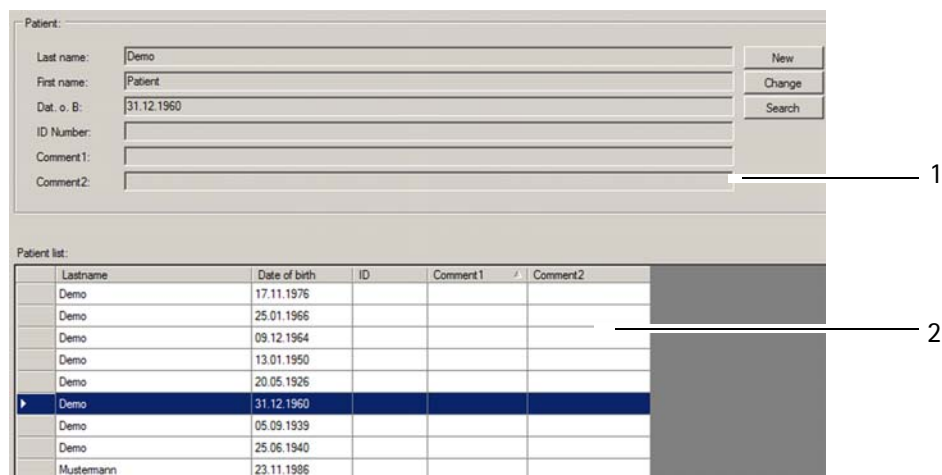

1 Attribute "Comment2" as an input field


2 Attribute "Comment2" in the patient list

Fig. 12-7: Custom attribute "Comment"



The deselection of attributes in the context menu is done in the same manner. The currently selected attributes are marked with a check.

12.5.2 "Devices" tab

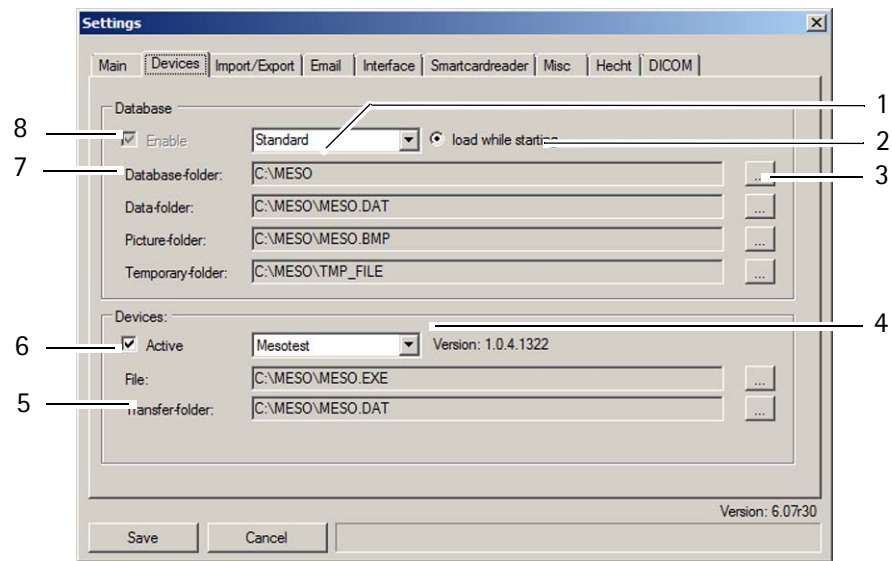

- | | |
|---|--|
| 1 Drop-down list for selecting a database | 6 Path for the currently selected device |
| 2 "Load while starting" radio button | 7 Checkbox for activating a device |
| 3 Folder selection | 8 Path for the database |
| 4 Drop-down list for selecting a device | 9 Checkbox for activating a database |

Fig. 12-8: "Devices" tab

"Database" box

Different users can be set up for different databases.

- ➔ In the drop-down list for selecting a database (1), choose the entry (User) that you would like to edit.
- ➔ Choose the corresponding path for the database, data, and images using the individual buttons for folder selection (3).
Normally, two different folders are created for data and images during installation, as shown here (7).
 - **For data:** Name of the device plus the .DAT identifier
 - **For images:** Name of the device plus the .BMP identifier
- ➔ For each user, enable whether the associated database should be enabled or not (8).

If more than one database is enabled, an additional drop-down list appears on the main page of the Patient Data Management user interface. You can enable a user (or the assigned database) in this list. The patient list and the associated examinations are updated when the active user is changed.

- ➔ You can also select the "load while starting" option (2) for exactly one user. The associated database is loaded by default when the Patient Data Management user interface is started, and is selected accordingly.

"Devices" box

Make settings for the connected devices in this box.

- ➔ In the drop-down list, select the desired device (4).
- ➔ If the device is in fact connected, then mark the checkbox [Active] (6).
- ➔ Choose the path to the associated device's application file using the folder selection button.


When you save patient and examination data in the database, it is initially placed in the so-called "Transfer folder". This folder is always created locally on the computer.

- ➔ Choose the transfer folder using the folder selection button. This should correspond with the name of the device and the .DAT extension.

You can also specify settings for devices which are not connected to the computer.

- ➔ Enable the checkbox [Write permission] for the devices that are in fact connected to the computer.

12.5.3 "Import/Export" tab

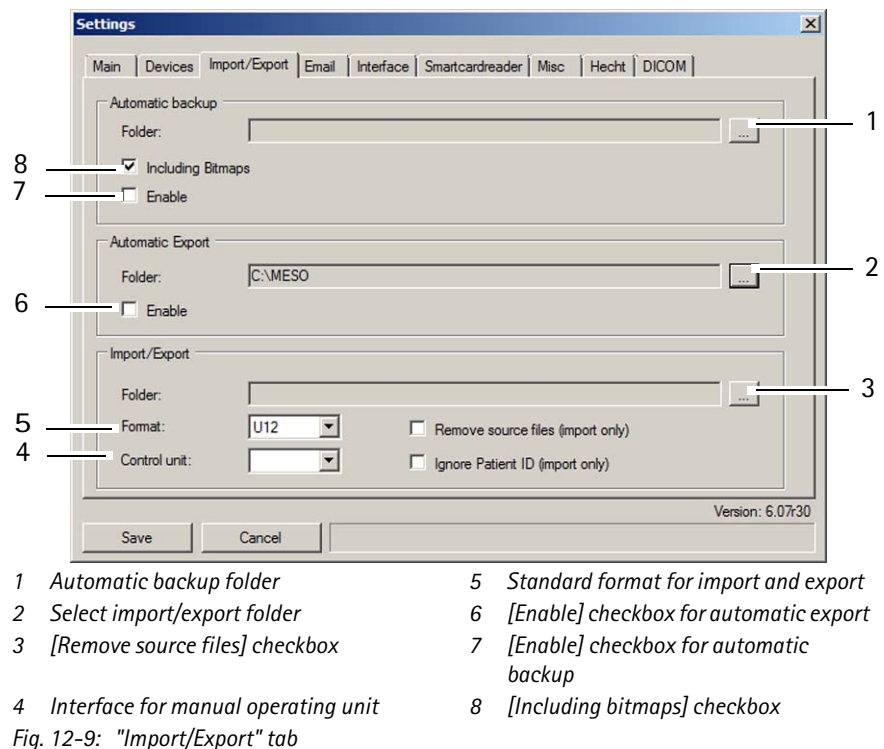

"Automatic backup" box

In addition to the manually executed backup ([sec. 12.4, page 55](#)), there is also the option to perform a backup when the Patient Data Management user interface is closed. The settings required for this can be specified in the group box.

- ➔ Select the folder where the data should be saved during an automatic backup using the folder selection button (1).
- ➔ Select the checkbox [Including Bitmaps] (8) if camera images should also be backed up.
- ➔ If the automatic backup should be executed with the specified settings, enable the [Enable] checkbox (7).

"Automatic export" box

If you enable this checkbox, every examination will be exported automatically.

- ➔ Select the folder using the folder selection button (2) which should normally be preselected for the automatic export.
- ➔ If the automatic export should be executed, enable the [Enable] checkbox (6).

"Import/Export" group box

You can enter settings in this box to import and export Patient Data Management data.




Note

The settings specified in this tab for import and export of data can be overwritten. These are only default values.

- ➔ Select the folder using the folder selection button (1) which should normally be used as the default for import or export.
- ➔ In the "Format" drop-down list (5), select whether the default import or export should be done for a folder or a single file (U12).

The checkbox „Control unit“ (4) is currently out of function.

12.5.4 "Email" tab

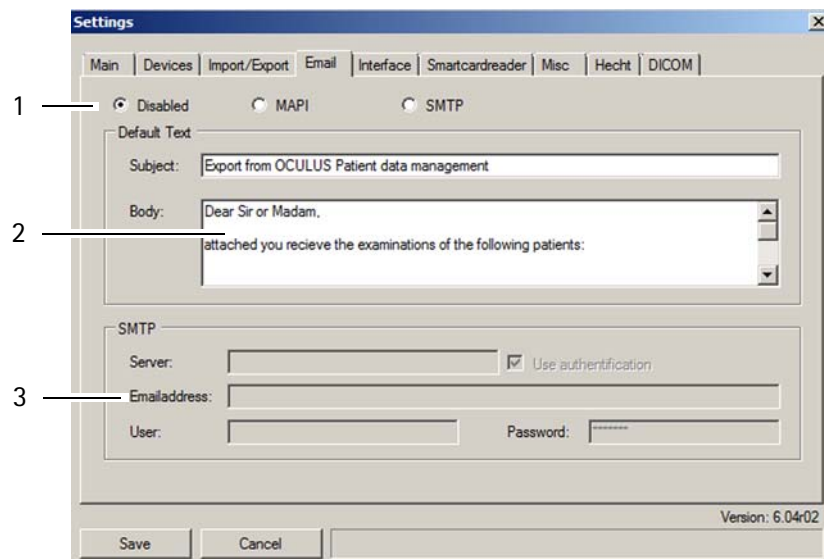

- 1 Select e-mail connection 3 SMTP connection selection
2 Standard text for subject and body text

Fig. 12-10: "Email" tab

With the three radio buttons in the upper part of the tab (1), you can determine whether the email connection is active, and if it is, how the transfer of data should take place.

- ➔ Enable the option "MAPI" if an email program (for example Microsoft Outlook) is installed on your computer and if data should be sent using this program.
- ➔ Enable the option "SMTP" if **no** email program is installed on your computer but you would still like to send the data by email. In this case, additional entries in the "SMTP" box are necessary.

"Default text" group box

- ➔ Enter the default text for sending emails in the "Subject" and "Body" (2) boxes. You can adjust this text before email is actually sent (for example particular to patients or examinations).

"SMTP" group box

If no email program is installed on the computer, several additional entries must be made here to be able to send emails.

- ➔ Contact your system administrator if you have questions about the individual entries.

12.5.5 "Interface" tab

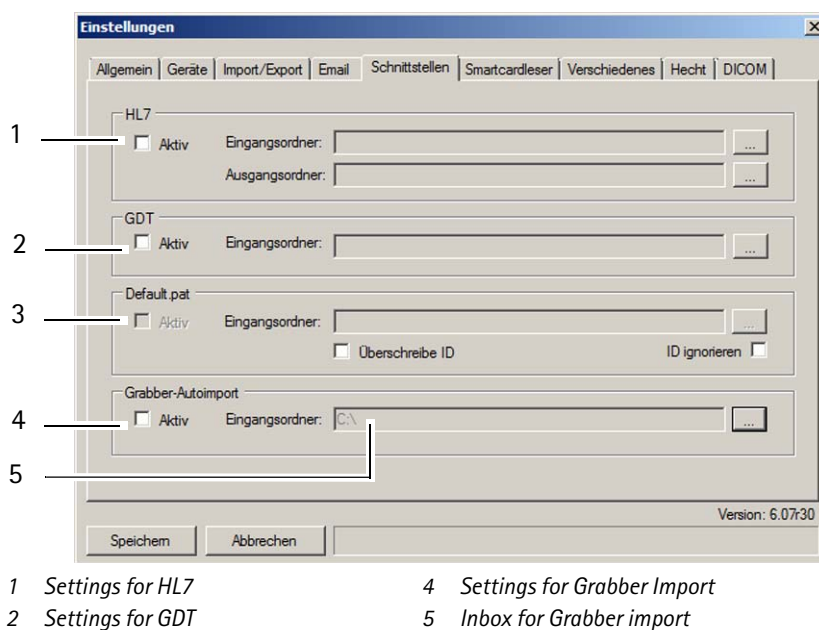


Fig. 12-11: "Interface" tab

Use this tab to enter the folder for different interface types and to enable or disable them.

The interface "Default.pat" (3) is used if the Patient Data Management user interface is launched by a third-party program.

- ➔ Add the interface "Default.pat" to the directory of the other program.
- ➔ Select the third-party program's folder as the "in" folder.

12.5.6 "Smartcardreader" tab


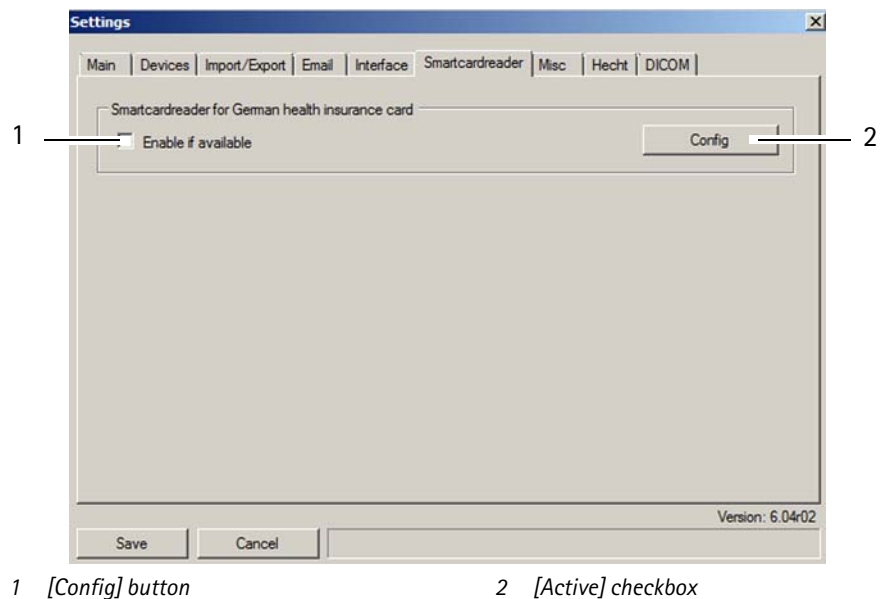


Fig. 12-12: "Smartcardreader" tab

You can configure a smartcard reader on this tab to be able to directly import patient data from the patient's insurance card into the Patient Data Management user interface. First, you have to configure the smartcard reader (usually only once).

- ➔ Press the [Config] button (2). A page appears for you to select the type of smartcard reader you are using.
- ➔ Select the [Enable if available] checkbox (1) to enable the smartcard reader.

12.5.7 "Misc" tab

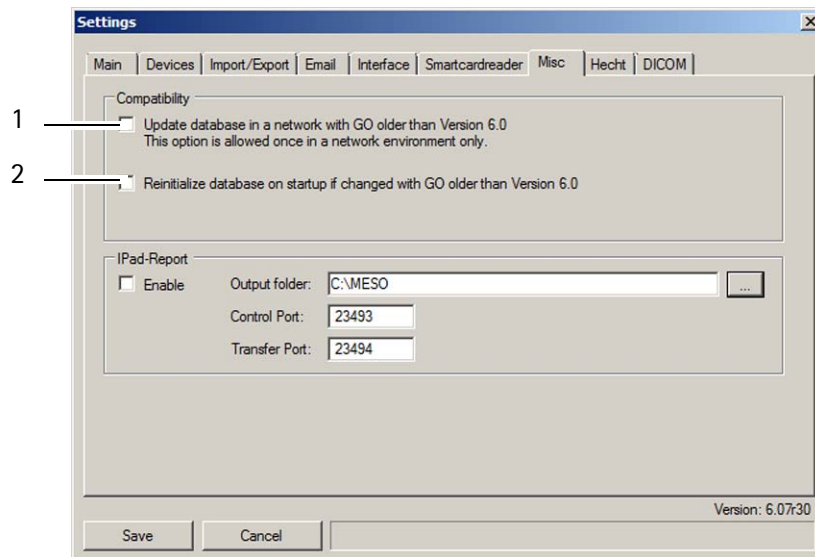

1 Checkbox to manage the database 2 Check box to import the database
Fig. 12-13: "Misc" tab

The Patient Data Management user interface is available in two different versions: V2.x and V6.x. In principle, implementing both versions on the same network should be avoided. In any case, if it is unavoidable due to the higher technical requirements of the new version V6.x, then the corresponding settings must be set on this tab.

- ➔ Enable the checkbox on this tab (1) on **exactly one** computer that has version V6.x installed.
- ➔ Make sure that this checkbox is **not** enabled on any other computer which has version V6.x of the Patient Data Management user interface installed.
- ➔ Enable Checkbox (2) so that data from the database will be read in again after the next start.

This ensures that the database is automatically updated when a patient is created on a computer with the older Patient Data Management V2.x. The "IPad Report" field is currently disabled.

12.5.8 "Hecht" and "DICOM" tabs

"Hecht" and "DICOM" tabs are currently disabled.

13 Cleaning, Disinfection and Maintenance

This chapter describes how to clean and disinfect the Mesotest II and how to change the viewer enclosure, the fuses and the panel illumination lamp.

To ensure satisfactory and reliable operation, we recommend that you have the Mesotest II checked every two years by our service department or an authorized dealer. If an error occurs which you are unable to correct, label the Mesotest II as "out of order" and contact our service department or an authorized dealer.



Attention


Risk of electric shock if the Mesotest II is not completely disconnected from the mains for these jobs.

- ➔ Turn the Mesotest II off, [sec. 7.5, page 19](#).
- ➔ Unplug the power cord before cleaning, disinfection and maintenance, [sec. 13.1, page 68](#).

- ➔ Always heed the product descriptions and directions for use of products you use to care for, clean, and disinfect the unit and/or its accessories.
- ➔ Do not clean the Mesotest II with aggressive, chlorinated, abrasive or sharp cleaning agents.

13.1 Unplug the Equipment

- ➔ Turn off the unit with the On/Off (3) switch.

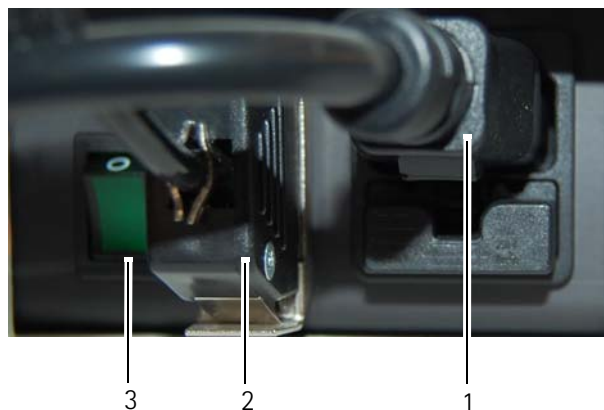

1 Mains connection

2 Plug for the control pad/computer/laptop

3 On/Off switch

Fig. 13-1: Plug connections of the Mesotest II

- ➔ Unplug the power cable from the socket (1).
- ➔ Push the button at the control pad plug and pull it out of the jack (2), or unplug the connecting cable to the computer/laptop.

13.2 Cleaning

Required Materials

- Cleaner for plastic surfaces with anti-static effect
- Cleaner for painted surfaces: Mixture of equal parts of alcohol and distilled water, possibly with a few drops of commercial detergent
- Soft cloth or lens brush
- Alcohol or lens cleaner
- Mild soap solution

Cleaning intervals: Wipe the viewer enclosure after each examination, and the housing whenever necessary.

- ➔ Turn the Mesotest II off, [sec. 7.5, page 19](#).
- ➔ Unplug the power cord.
- ➔ When cleaning, use a damp cloth and make sure that no liquid enters the Mesotest II.
- ➔ Clean the plastic surfaces and painted surfaces with appropriate cleaning agents.
- ➔ Clean the lenses with a soft cloth or lens brush, and, if necessary, with alcohol or a lens cleaner.
- ➔ Clean the viewer enclosure with a mild soap solution.

13.3 Disinfection

Required Materials

To disinfect all surfaces (except the Plexiglas back pane of the viewer) we recommend the use of:

- Mikrozid sensitive wipes premium
Make: Schülke & Mayr GmbH; Softpack 48x, Art.No. 165711
- ➔ Disinfect the viewer enclosure after every examination and the housing whenever necessary.

If the viewer enclosure is very dirty, you can exchange it, [sec. 13.4.1, page 71](#).



Note:

Equipment damage caused by disinfectant solution

The disinfectant solution may damage the finish if it is sprayed directly on it.

- ➔ Spray the disinfectant solution onto a cleaning cloth, do not spray it directly on the device.

13.4 Care and Maintenance

To ensure satisfactory and reliable operation, we recommend that you

- Have the Mesotest II checked every two years by our service department or an authorized dealer.



Note:

Faulty examinations due to damaged equipment

If an error occurs that you can not resolve

- Label the damaged Mesotest II as "out of order".
 - Report the damage to OCULUS Service or your authorized dealer.
 - Use only undamaged Mesotest II devices.
-

13.4.1 Exchange Viewer Enclosure

If the viewer enclosure (1) is worn or is very dirty, you can exchange it.

Required Materials

- Soft cloth as underlay
- New viewer enclosure (Order number 016280001003)

Proceed as follows:

- ➔ Turn the Mesotest II off, [sec. 7.5, page 19](#).
- ➔ Unplug the power cord.
- ➔ Lay the Mesotest II on its back.



Note:

Risk of damaging the viewer glass if the unit is put down incorrectly

- ➔ Lay the Mesotest II on its back only. Use the soft cloth as an underlay.

Removal:

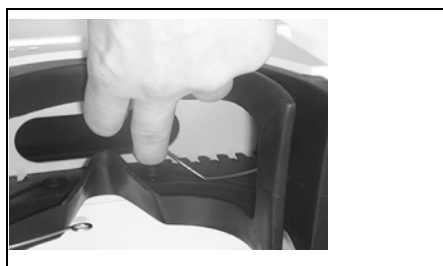


Fig. 13-2: Removing the yoke springs

- ➔ Remove the yoke springs:
 - To do this, take hold of the yoke springs one after the other at the two round ends and unhook the spring ends from the rubber part by pulling them outwards.
 - Then fully remove the yoke spring from the top guide groove.



Fig. 13-3: Removing the viewer enclosure

- ➔ Then press the viewer enclosure out of its guide groove.
- ➔ Take hold of one side of the viewer enclosure and pull it out.

Insertion:

- ➔ Compress the viewer enclosure vertically and insert it into the frame of the housing halfshell.
- ➔ Starting at one side, press the groove of the viewer enclosure into the frame.

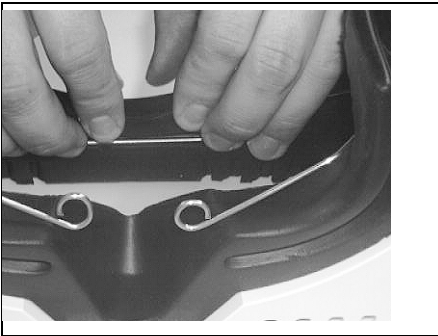


Fig. 13-4: Inserting the yoke springs

- Make sure that the small air ducts on the back are not folded over.
- First insert continuous yoke of the yoke spring into the top guide groove.
- Make sure that the yoke sits fully in the groove.
- Then take hold of the two bent ends of the yoke spring, one after the other, and hook them back into place in the bottom groove.

13.4.2 Replacing the Fuse

The Mesotest II has two fuses. These are located in a small fuse box integrated in the mains power plug.

You can replace a blown fuse.



Note:

Function loss due to incorrect fuse

- Use only the fuse type that is specified on the rating plate, [sec. 2, page 2](#).

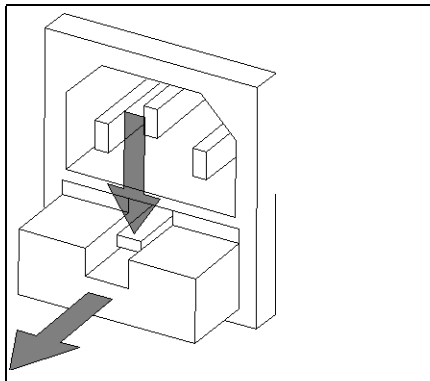


Fig. 13-5: Pull out the fuse holder

- Turn the Mesotest II off, [sec. 7.5, page 19](#).
- Unplug the power cord.
- Open the drawer. To do this press on the lug on the top of the drawer.
- Pull out the drawer.
- Replace the blown fuse.
- You can recognise a blown fuse by the burnt filament.
- Push the drawer back in until the lug snaps back into place.
- Connect the Mesotest II to the mains.
- You can now turn on the Mesotest II and start doing exams.

13.4.3 Change the Lamp at the Panel Illumination

You can change a faulty lamp.



Note:

Risk of functional damage caused by wrong lamp.

→ Use only genuine lamps from Oculus. Order number 016280006004.

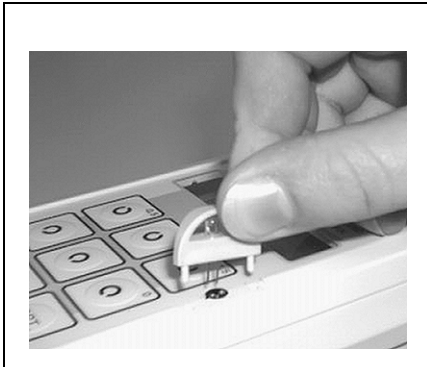


Fig. 13-6: Changing the lamp

- The bulb is changed together with the plastic holder, which also acts as a glare shield.
- Simply pull the lamp up and remove it.
- Insert the lamp.
When inserting the new lamp, take care not to bend the pins (connecting wires) of the bulb.
- Push the lamp into its fixture.

14 Troubleshooting



Attention

Risk of personal injury or equipment damage due to improper trouble-shooting

- If an error occurs which you are unable to correct by following the instructions below, label the device as „out-of-order“ and contact our service department or an authorized dealer.

Error	Possible source	Remedy
Nothing happens when the power switch is turned on	The Mesotest II is not connected to the power supply	Plug the power cable into the socket or into the low heat device socket at the Mesotest II.
	Power failure or power socket is not active	Inform the in-house electrician
	The control pad or the serial cable of the computer is not connected properly	Check that the connector is plugged in properly
	The device fuses have blown	Change the fuses (see sec. 13.4.2, page 72)
	The unit has been switched off and back on again too quickly	Wait approx. 5 seconds before turning the unit back on again
Patient claims he cannot see anything	The unit is in stand-by mode	Press any button on the control pad
	The examination room is too bright	Darken the room or cover the back pane of the viewer with the supplied cover.
Panel illumination is not functioning	Lamp at the panel illumination is faulty	Change the lamp (see sec. 13.4.3, page 73)
R/L/BIN and -DPT. with corresponding minus lens powers are not working despite display in the LCD	Supplementary equipment is not installed in the device	

15 Transport and Storage

The Mesotest II, must be properly dismantled and packed before being transported or stored.

15.1 Transport and Storage Information

Storage

This device does not keep up with the temperature conditions for storage according to ISO 15004-1.

- ➔ Do not store the Mesotest II under conditions over 50 °C and under -10 °C.

In accordance with DIN EN 60601-1, the storage conditions are:

Ambient temperature range	-10°C to +55°C
Relative humidity, including condensation	10% to 95%
Air pressure range	700 hPa to 1060 hPa

Transport

In accordance with DIN EN 60601-1, the transport conditions are:

Ambient temperature range	-40°C to +70°C
Relative humidity, including condensation	10% to 95%
Air pressure range	500 hPa to 1060 hPa

After storage and transport

- ➔ Wait approx. 3-4 hours after transport before operating the Mesotest II for the first time. In the event of extreme temperature changes from cold areas to warm rooms, the optical components can become fogged.

15.2 Disassembly and Packing

- End the current session.
- Pull the power plug and disconnect the plug connectors.



Attention

Equipment damage due to improper transport and storage


- Avoid knocks, vibrations, dirt, high temperatures and moisture.

-
- Transport the Mesotest II carefully.
 - Avoid shocks, vibration and dust.
 - Store the Mesotest II in accordance with the storage conditions.
Avoid storing in the vicinity of heaters and moisture.



If you want to send the Mesotest II to OCULUS for repair, you can order an empty packaging box from OCULUS Service. That will help prevent potential transport damage.

16 Disposal

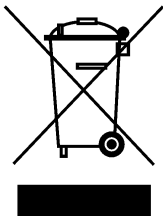

In accordance with Guideline 2002/96/EG of the European Parliament and of the Council dated 27 January 2003, and also the Law of the Federal Republic of Germany on the Commercialization, Recall and Environmentally Compatible Disposal of Electrical and Electronic Equipment, old electrical and electronic equipment must be sent out for recycling and may not be disposed in household trash.

- Transport the Mesotest II in a compliant manner.

17 Terms of Warranty and Service

Any software included in the delivery was tested by us and complies with technical standards. Note the following warranty provisions:

- Prior to and while operating the device it is important that you observe the instruction manual and safety instructions.
- The Mesotest II carries a warranty to which you are entitled in accordance with the legal provisions.
- If any unauthorized persons interfere with the Mesotest II, all warranty entitlements shall be void. Any inappropriate modifications or repairs can cause grave danger to the user and patient.
- Any entitlement to a warranty shall also be void if unauthorized persons interfere with the computer hardware and supplied software.
- Make transport damage claims to the shipping company during or immediately after delivery. Have the damage confirmed on the bill of lading, so that a proper claim settlement is possible.
- In general, the general terms and conditions of business and delivery apply as per the date of purchase.

17.1 Liability for proper function or damages

OCULUS will only accept responsibility for the safety, reliability and serviceability of the Mesotest II if the unit is used in compliance with the following terms:

- Use the device in accordance with these instructions and the accompanying user manual.
- There are no user-serviceable parts either on or inside the Mesotest II. OCULUS shall not assume any liability if assembly, extensions, adjustments, changes or repairs are carried out by unauthorized personnel; if the Mesotest II is maintained improperly; or if it is handled incorrectly.
- If the work described above is carried out by persons authorized to do so, they must be required to supply a certification detailing the nature and scope of repairs, and, if applicable, to specify modifications to the rated data and area of work. The certificate must bear a date, a signature, specify who carried out the work, and contain company information.
- On request, and for this purpose, OCULUS will supply authorized persons with spare parts lists and additional descriptions.
- Make sure that only original OCULUS parts are used for service and maintenance.

17.2 Address of the manufacturer and service department

Our service department or authorized representatives will furnish you with additional information. Address of the manufacturer and service department:

Germany:

OCULUS Optikgeräte GmbH
Münchholzhäuser Straße 29
35582 Wetzlar, Germany
Tel.: 06412005-0
Fax: 06412005-255
E-mail: sales@oculus.de
www.oculus.de



USA:

OCULUS, Inc.
17721 59th Avenue NE
Arlington
WA 98223-1337
Tel. +1 425-670-9977
Fax +1 425-670-0742
e-mail: sales@oculususa.com
<http://www.oculususa.com>



18 Declaration of Conformity



OCULUS Optikgeräte GmbH
Münchholzhäuser Str. 29
D-35582 Wetzlar
Germany

Tel: +49 (0)641 / 20 05 – 0
Fax: +49 (0)641 / 20 05 – 255

Konformitätserklärung Declaration of Conformity

Wir, OCULUS Optikgeräte GmbH, erklären in alleiniger Verantwortung, dass die unten beschriebenen Medizinprodukte allen Anforderungen der nachstehenden Richtlinie und den damit verbundenen harmonisierten Normen entsprechen: 93/42/EWG

We, OCULUS Optikgeräte GmbH, declare on our own responsibility that the medical devices described below are in compliance with requirements of the following directive and their related harmonized standards: 93/42/EEC

**Produktbezeichnung/
Product name**

Sehtestgerät/
Visual acuity tester

**Artikelnummer und Typ/
Article number and type**

62800 Mesotest II
62801 Mesotest II konkav

**Konformitätsbewertungsverfahren/
Conformity assessment procedure**

Richtlinie 93/42/EWG: Anhang VII
Conformity according: 93/42/EEC, Annex VII

**Ort, Datum/
Place, date**

Wetzlar, 26.02.2014

**Name und Funktion/
Name and function**

Qualitätsmanagement / Quality Manager
OCULUS Optikgeräte GmbH



Eckhard Loh

19 Technical Data

Technical specifications

Dimensions (W x D x H)	246 x 377 x 464 mm (9.7 x 14.9 x 18.2 inches)
Weight (without control pad and power cable)	6.5 kg1 (4.3 lbs)
Max. power consumption	36 VA
Voltage	90-264 V AC
Frequency	50 / 60 Hz
Control pad:	
Dimensions (W x D x H)	82 x 152 x 45 mm (3.2 x 6.0 x 1.7 inches)
Weight	256 gr (0.6 lbs)
Shelf life	up to 10 years

Ambient conditions

Temperature	+10°C to +35°C
Humidity	30% to 90%
Air pressure	800 hPa to 1060 hPa

Storage conditions according to DIN EN 60601-1

Ambient temperature range	-10°C to +55°C
Relative humidity, including condensation	10% to 95%
Air pressure range	700 hPa to 1060 hPa

Transport conditions according to DIN EN 60601-1

Ambient temperature range	-40°C to +70°C
Relative humidity, including condensation	10% to 95%
Air pressure range	500 hPa to 1060 hPa

Classification according to DIN EN 60601-1

Type of protection against electric shock	Protection class 1
Level of protection against electric shock	Type B
Level of protection against harmful penetration of water	IP 20

Computer

Use a computer which is in conformity with the DIN EN 60950 standard.

Recommended computer specifications (for optional software)	2 GB RAM, 320 GB HDD, 1.8 GHz, Windows® 8, 64bit
---	--

20 Appendix

20.1 Electromagnetic Compatibility

Medical electrical equipment is subject to special precautionary requirements with respect to EMC, and must be installed and operated according to the EMC-Instructions contained in the accompanying paperwork.

No special measures need be observed in respect of OCULUS devices and systems.

Portable and mobile RF-communications devices can interfere with electrically operated medical devices.

Minimal performance quality and essential performance criteria

- A short flicker of the illumination of the Mesotest II during the examination is permissible because it will not affect the functionality described in the intended use.
- A short interruption of the USB connection during the examination is permissible because it will not affect the functionality described in the intended use.



Attention

The use of accessories, transducers and cables not specified by OCULUS (for example as replacement parts) may result in increased emissions or decreased immunity of the Mesotest II.

- ➔ Use only the original accessories, transducers and cables specified by OCULUS.

The use of accessories, transducers and cables specified by OCULUS with devices other than the Mesotest II may result in increased emissions or decreased immunity of the other device.

- ➔ Do not use the accessories, transducers and cables specified by OCULUS with devices other than the Mesotest II.
-

To be in compliance with the requirements of the IEC 60601-1-2: 6.1 and 6.2 the following types of equipment, accessories, power adapters and cables must be used.

Order number	Description	
62800	Mesotest II	
62801	Mesotest II konkav	
05200320	Kabel with plug, EU standard	2.5 m (8.2021 feet)
05200210 (110 Volt)	Kabel with plug, US standard	2.5 m (8.2021 feet)

20.2 Guidance and Manufacturer's Declaration - Electromagnetic Emissions for the Mesotest II

Guidance and manufacturer's declaration electromagnetic emissions
IEC 60601-1-2, 5.2.2.1, table 1

The OCULUS Mesotest II is intended for operation in the electromagnetic environment specified below. The user of the Mesotest II should ensure that it is being used in such an environment.

Emissions test	Compliance	Electromagnetic environment – guidance
RF emissions CISPR 11	Group 1	The Mesotest II uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
HF-emissions CISPR 11	Class B	
Harmonics emissions IEC 61000-3-2	Class A	
Voltage fluctuations / flicker emissions IEC 61000-3-3	complies	

Guidance and manufacturer's declaration electromagnetic immunity,
IEC 60601-1-2, 5.2.2.1, table 2

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
Electrostatic discharge (ESD) IEC 61000-4-2	± 6 kV contact ± 8 kV air	± 6 kV contact ± 8 kV air	Floors should be made of wood or concrete or covered with ceramic tiles. If the floor is covered with synthetic material, the relative humidity must be at least 30%.
Electrical Fast transient/bursts IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV for input/output lines	± 2 kV for power supply lines ± 1 kV for input/output lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 6100-4-5	± 1 kV line(s) to line(s) ± 2 kV line(s) to earth	± 1 kV differential mode ± 2 kV line(s) to earth	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions, and voltage variations on power supply input lines IEC 61000-4-11	< 5% U_T (> 95% dip in U_T) for 0,5 cycle 40 % U_T (60% dip in U_T) for 5 cycles 70% U_T (30% dip in U_T) for 25 cycles <5% U_T (> 95% dip in U_T) for 5 s	< 5% U_T (> 95% dip in U_T) for 0,5 cycle 40 % U_T (60% dip in U_T) for 5 cycles 70% U_T (30% dip in U_T) for 25 cycles <5% U_T (> 95% dip in U_T) for 5 sec	Mains power quality should be that of a typical commercial or hospital environment. If the user of the Mesotest II requires continued operation during power mains interruptions, it is recommended that the Mesotest II be powered from an uninterruptible power supply or battery.
Power frequency (50/60Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
Note: U_T is the a.c. mains voltage prior to application of the test level.			

Guidance and manufacturer's declaration electromagnetic immunity,
IEC 60601-1-2, 5.2.2.2, table 4

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment – Guidelines
			Portable and mobile RF communications equipment should be used no closer to any part of Mesotest II, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance
Conducted RF IEC 61000-4-6	3 Vrms 150 KHz to 80 Mhz	$V_{rms} = 3 \text{ V}$	$d = \left[\frac{3,5}{(V_1)} \right] \sqrt{P}$
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.5 GHz	$E = 3 \text{ V/m}$	$d = \left[\frac{3,5}{(E_1)} \right] \sqrt{P} \quad 80\text{MHz to } 800 \text{ MHz}$ $d = \left[\frac{7}{(E_1)} \right] \sqrt{P} \quad 800 \text{ MHz to } 2.5 \text{ GHz}$ where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m). Field strength from fixed RF transmitters, as determined by an electromagnetic site survey (a), should be less than the compliance level in each frequency range (b). Interface may occur in the vicinity of equipment marked with the following symbol: 
Note 1:	At 80 Hz and 800 MHz, the higher frequency range applies.		
Note 2:	These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.		
a. Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radios, AM and FM radio broadcast, and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, and electromagnetic site survey should be considered. If the measured field strength in the location in which the Mesotest II is used exceeds the applicable RF compliance level above, the Mesotest II should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the Mesotest II.			
b. Over the frequency range 150 KHz to 80 MHz, field strengths should be less than 3 V/m.			

Recommended separation distances between portable and mobile RF communications equipment and the Mesotest II, IEC 60601-1-2, 5.2.2.2, table 6

The Mesotest II is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the Mesotest II can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the Mesotest II as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m		
	150 KHz to 80 Mhz $d = 1.2 \sqrt{P}$	80 MHz to 800 MHz $d = 1.2 \sqrt{P}$	800 MHz to 2.5 GHz $d = 2.3 \sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.80	3.80	7.3
100	12	12	23

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in metres (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.