


OCULUS Trial Lens Cases for t1s Trial Lenses Modular Concept for Professional Refraction




OCULUS Trial Lenses

Clever detail solutions that impress

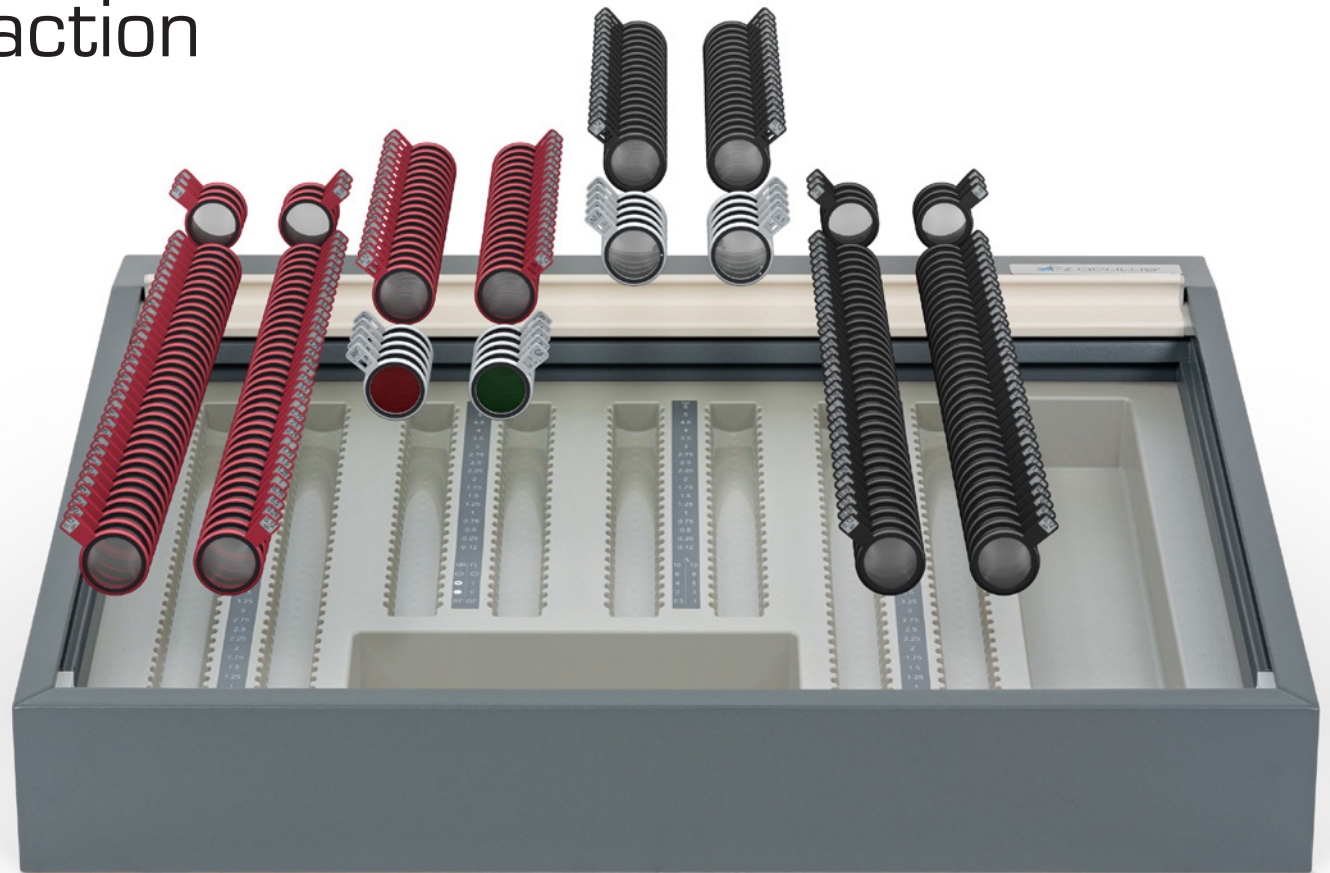

Tab with aperture and reading mark for precise alignment

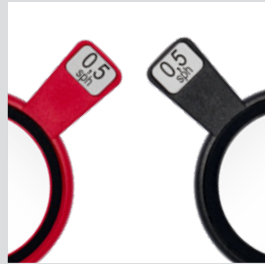

Clear colour coding of concave and convex trial lenses




Clear indication of the axis position of cylindrical trial lenses



Clear indication of the base position of prism trial lenses

Important information

The tls trial lenses* (tls = trial lens small) do have a rim diameter of 28 mm (lens diameter 21 mm).

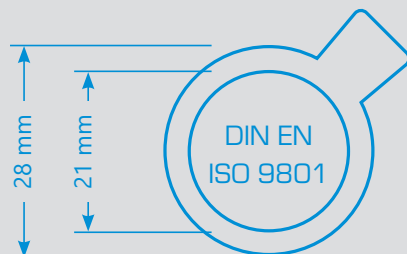

Coating of the trial lenses (optional)

Super-broadband antireflection coating with residual green reflection and hard coating for optimal resistance and longevity.


No annoying "snagging"

All OCULUS tls trial lenses are plano-convex or plano-concave grounded, enabling the trial lenses to be switched out effortlessly.

Exclusively at OCULUS!

All tls trial lenses can be ordered individually from your local OCULUS representative. You can find them at

www.oculus.de/en/distributor-locator

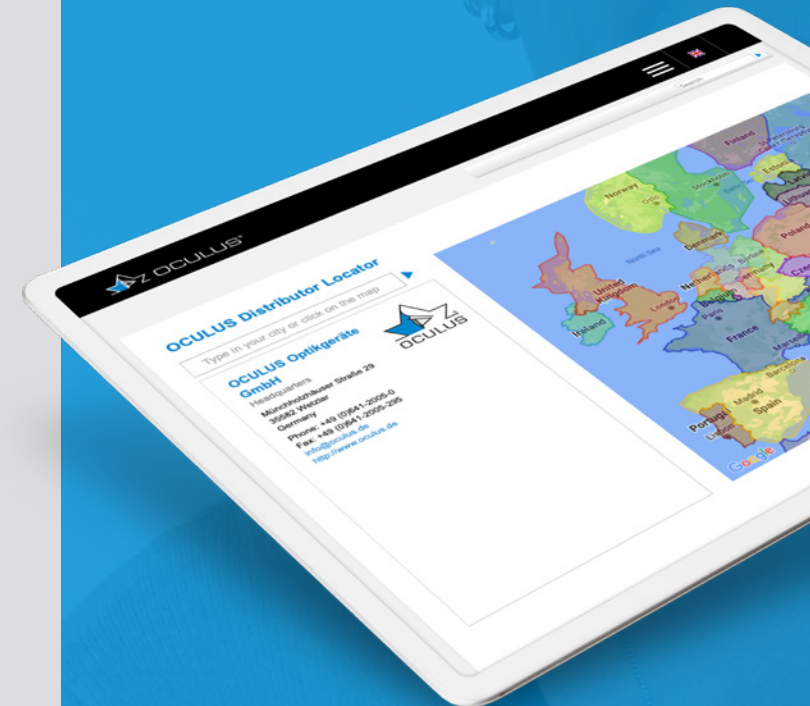

* They are suited for the OCULUS Universal Trial Frame UB5.



OCULUS Trial Lens Cases with Modular Concept

Unrivalled quality made by OCULUS. Flexible, sustainable and individually expandable to ensure lively variety.

Your new trial lens case is configured the way you want and meets your own personal requirements:

[Preconfigured trial lens sets and options](#)

Rethought and redesigned. User-friendly solutions with growth potential.

[Modules](#)

Use the modules to set up your personalized refraction equipment.

[Single trial lenses](#)



By enabling you to order individual lenses we offer you greater flexibility and more sustainability.

[Trial lens cases and lens trays](#)

Sophisticated ergonomic storage solutions round our concept.

[We inspire, you define!](#)

The new OCULUS concept is based on the motto "Only buy what you truly need".

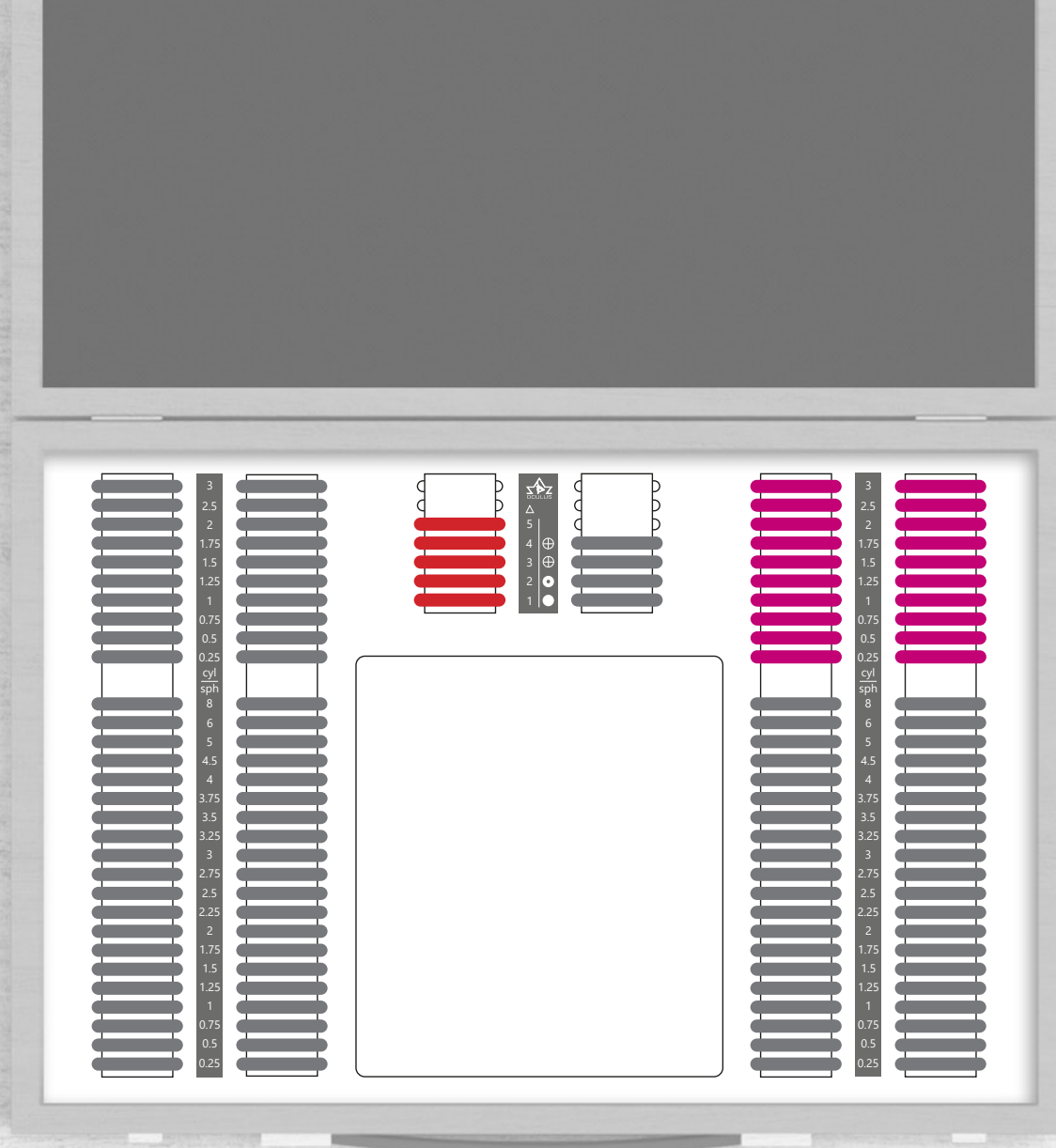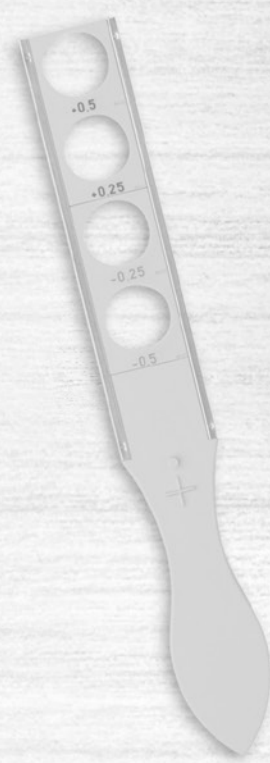

378 mm (14.9 in)

255 mm (10.0 in)



Trial Lens Set Basic S

The starter

The following modules are already included in the Trial Lens Set Basic S

Description	Article no.	Quantity
tls trial lenses +/- sph (basic)	10033238	80
tls trial lenses - cyl (basic) <i>axis position 90°</i>	10033247	20
tls trial lenses special lenses (basic) <i>2 x cross hair, 1 x occluder, 1 x pinhole disc 1.25</i>	10033256	4
Trial Lens Set Basic S*	10033416	104

* uncoated



We recommend the following optional modules for expansion

Description	Article no.	Quantity
tls trial lenses + cyl (basic) <i>axis position 90°</i>	10033249	20
tls trial lenses prism. (basic) <i>base 0°</i>	10033259	5

Everything perfectly organized!

You will find the ideal storage solution for your new trial lenses on **page 20**.

The illustration above gives a schematic overview of the trial lens case (article no. 10035803) with the Trial Lens Set Basic S and the optional modules.

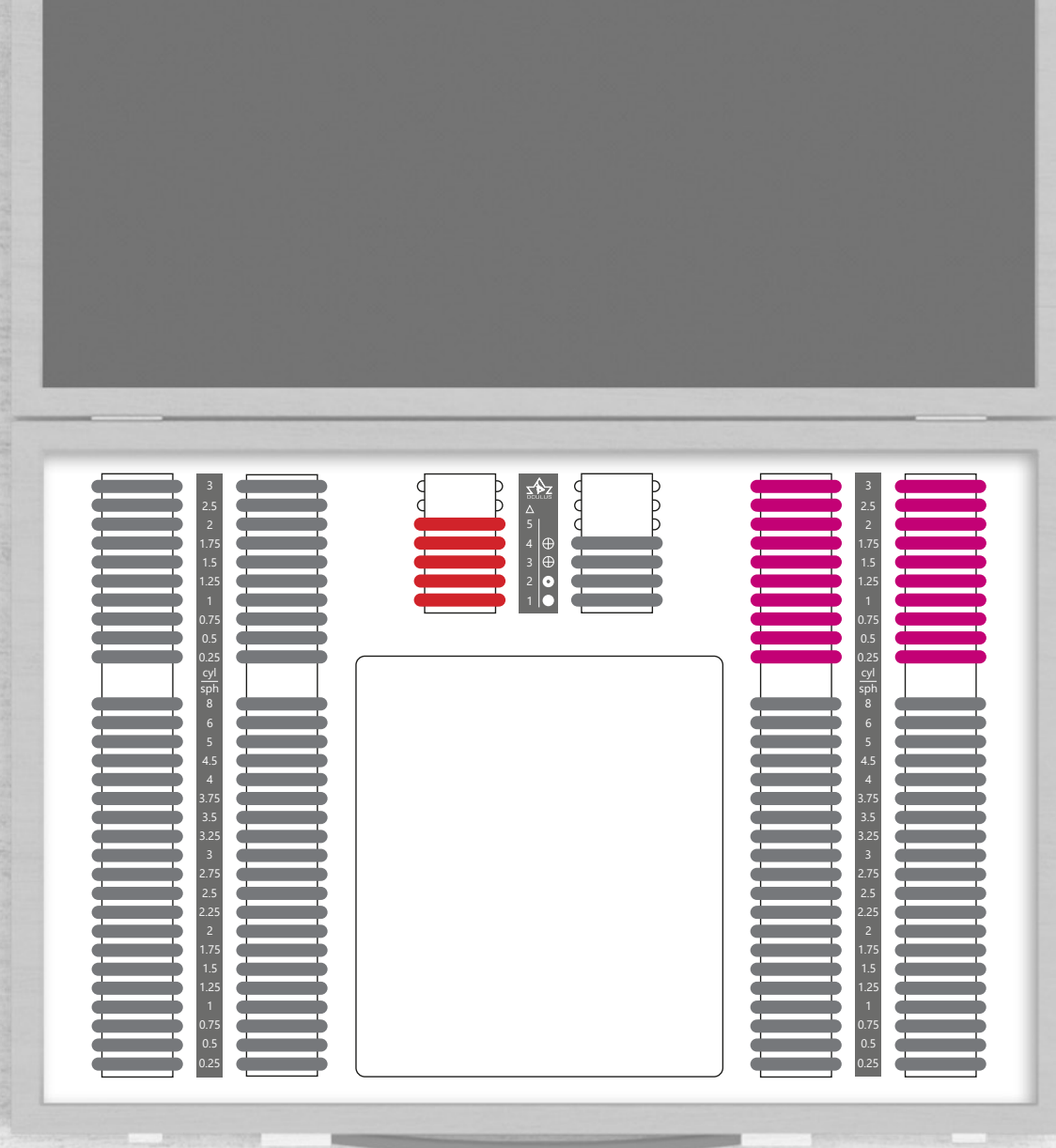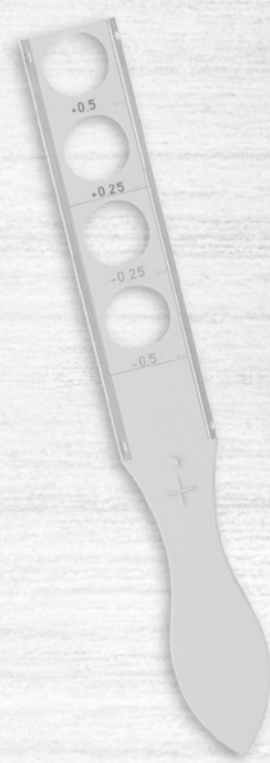

378 mm (14.9 in)

255 mm (10.0 in)



Trial Lens Set Basic S Coated

The starter

The following modules are already included in the Trial Lens Set Basic S coated

Description	Article no.	Quantity
tls trial lenses +/- sph (basic), coated	10033240	80
tls trial lenses - cyl (basic), coated <i>axis position 90°</i>	10033248	20
tls trial lenses special lenses (basic) <i>2 x cross hair, 1 x occluder, 1 x pinhole disc 1.25</i>	10033256	4
Trial Lens Set Basic S, coated	10033417	104







We recommend the following optional modules for expansion

Description	Article no.	Quantity
tls trial lenses + cyl (basic), coated <i>axis position 90°</i>	10033250	20
tls trial lenses prism. (basic), coated <i>base 0°</i>	10033260	5

Everything perfectly organized!

You will find the ideal storage solution for your new trial lenses on **page 20**.

The illustration above gives a schematic overview of the trial lens case (article no. 10035803) with the Trial Lens Set Basic S, coated and the optional modules.

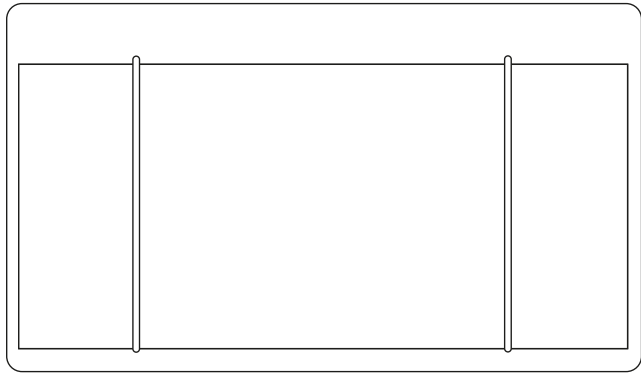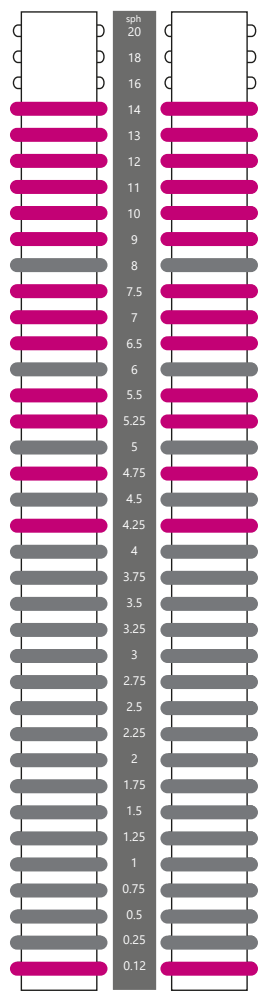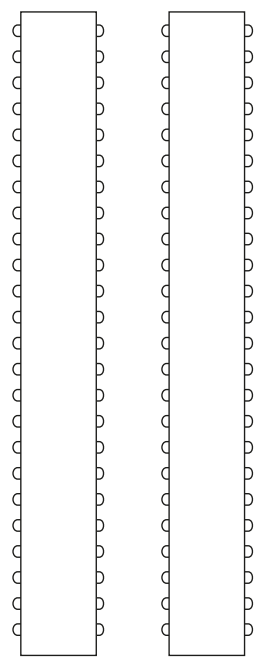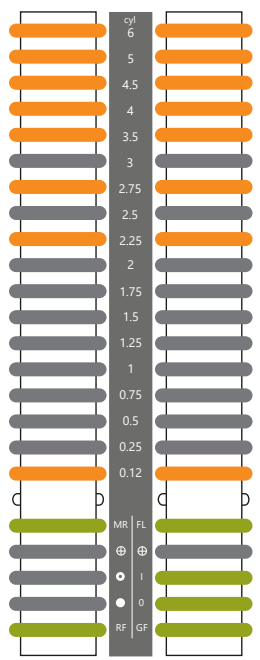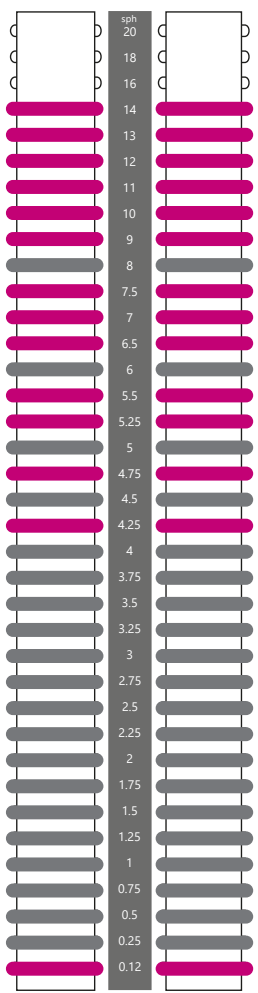

476 mm (18.7 in)

578 mm (22.8 in)

Trial Lens Set Basic S Coated

The starter in Trial Lens Case Comfort S and Premium S

The following modules are already included in the Trial Lens Set Basic S coated

Description	Article no.	Quantity
tls trial lenses +/- sph (basic), coated	10033240	80
tls trial lenses - cyl (basic), coated <i>axis position 90°</i>	10033248	20
tls trial lenses special lenses (basic) <i>2 x cross hair, 1 x occluder, 1 x pinhole disc 1.25</i>	10033256	4
Trial Lens Set Basic S, coated	10033417	104

We recommend the following optional modules for expansion

Description	Article no.	Quantity
tls trial lenses special lenses (extension) <i>coloured lens red, coloured lens green, frosted lens, plano lens, stenopaic slit 1.25 mm, Maddox cylinder red</i>	10033257	6
tls trial lenses +/- sph (extension), coated	10033242	56
tls trial lenses - cyl (extension), coated <i>axis position 90°</i>	10033252	16

Everything perfectly organized!

You will find the ideal storage solution for your new trial lenses on **page 20**.



The illustration above gives a schematic overview of the trial lens case (article no. 10035809) with the Trial Lens Set Basic S, coated and the optional modules.

20	sph	20
18		18
16		16
14		14
13		13
12		12
11		11
10		10
9		9
8		8
7.5		7.5
7		7
6.5		6.5
6		6
5.5		5.5
5.25		5.25
5		5
4.75		4.75
4.5		4.5
4.25		4.25
4		4
3.75		3.75
3.5		3.5
3.25		3.25
3		3
2.75		2.75
2.5		2.5
2.25		2.25
2		2
1.75		1.75
1.5		1.5
1.25		1.25
1		1
0.75		0.75
0.5		0.5
0.25		0.25
0.12		0.12

6	cyl	6
5		5
4.5		4.5
4		4
3.5		3.5
3		3
2.75		2.75
2.5		2.5
2.25		2.25
2		2
1.75		1.75
1.5		1.5
1.25		1.25
1		1
0.75		0.75
0.5		0.5
0.25		0.25
0.12		0.12
MR	FL	
⊕	⊕	
⊖	⊖	
0		
RF	GF	

6	cyl	6
5		5
4.5		4.5
4		4
3.5		3.5
3		3
2.75		2.75
2.5		2.5
2.25		2.25
2		2
1.75		1.75
1.5		1.5
1.25		1.25
1		1
0.75		0.75
0.5		0.5
0.25		0.25
0.12		0.12
Δ		
10		12
6		8
4		5
2		3
0.5		1

20	sph	20
18		18
16		16
14		14
13		13
12		12
11		11
10		10
9		9
8		8
7.5		7.5
7		7
6.5		6.5
6		6
5.5		5.5
5.25		5.25
5		5
4.75		4.75
4.5		4.5
4.25		4.25
4		4
3.75		3.75
3.5		3.5
3.25		3.25
3		3
2.75		2.75
2.5		2.5
2.25		2.25
2		2
1.75		1.75
1.5		1.5
1.25		1.25
1		1
0.75		0.75
0.5		0.5
0.25		0.25
0.12		0.12

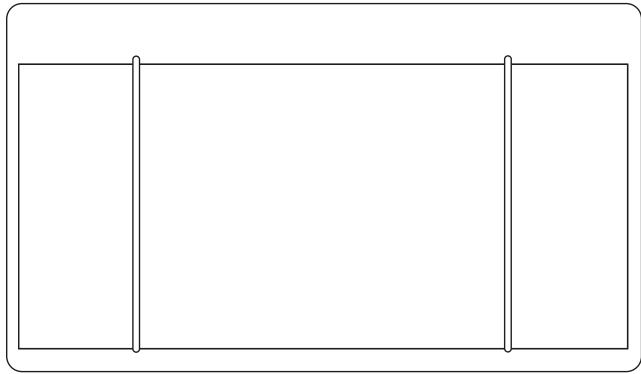

476 mm (18.7 in)

578 mm (22.8 in)

Trial Lens Set Comfort S

The classic

The following modules are already included in the Trial Lens Set Comfort S

Description	Article no.	Quantity
tls trial lenses +/- sph (basic)	10033238	80
tls trial lenses +/- sph (extension)	10033241	56
tls trial lenses - cyl (basic) <i>axis position 90°</i>	10033247	20
tls trial lenses - cyl (extension) <i>axis position 90°</i>	10033251	16
tls trial lenses special lenses (basic) <i>2 x cross hair, 1 x occluder, 1 x pinhole disc 1.25</i>	10033256	4
tls trial lenses special lenses (extension) <i>coloured lens red, coloured lens green, frosted lens, plano lens, stenopaic slit 1.25 mm, Maddox cylinder red</i>	10033257	6
Trial Lens Set Comfort S*	10033418	182

* uncoated

We recommend the following optional modules for expansion

Description	Article no.	Quantity
tls trial lenses +/- sph <i>each 2 x 16/18/20 D</i>	10033243	12
tls trial lenses + cyl (basic) <i>axis position 90°</i>	10033249	20
tls trial lenses + cyl (extension) <i>axis position 90°</i>	10033253	16
tls trial lenses prism. (basic) <i>base 0°</i>	10033259	5
tls trial lenses prism. (extension) <i>base 0°</i>	10033261	5

Everything perfectly organized!

You will find the ideal storage solution for your new trial lenses on **page 20**.



The illustration above gives a schematic overview of the trial lens case (article no. 10035809) with the Trial Lens Set Comfort S and the optional modules.

20	18	16	14	13	12	11	10	9	8	7.5	7	6.5	6	5.5	5.25	5	4.75	4.5	4.25	4	3.75	3.5	3.25	3	2.75	2.5	2.25	2	1.75	1.5	1.25	1	0.75	0.5	0.25	0.12
----	----	----	----	----	----	----	----	---	---	-----	---	-----	---	-----	------	---	------	-----	------	---	------	-----	------	---	------	-----	------	---	------	-----	------	---	------	-----	------	------

6	5	4.5	4	3.5	3	2.75	2.5	2.25	2	1.75	1.5	1.25	1	0.75	0.5	0.25	0.12	MR	FL	⊕	⊕	⊖	⊖	0	RF	GF
---	---	-----	---	-----	---	------	-----	------	---	------	-----	------	---	------	-----	------	------	----	----	---	---	---	---	---	----	----

6	5	4.5	4	3.5	3	2.75	2.5	2.25	2	1.75	1.5	1.25	1	0.75	0.5	0.25	0.12	Δ	12	10	8	6	5	4	3	2	1
---	---	-----	---	-----	---	------	-----	------	---	------	-----	------	---	------	-----	------	------	---	----	----	---	---	---	---	---	---	---

20	18	16	14	13	12	11	10	9	8	7.5	7	6.5	6	5.5	5.25	5	4.75	4.5	4.25	4	3.75	3.5	3.25	3	2.75	2.5	2.25	2	1.75	1.5	1.25	1	0.75	0.5	0.25	0.12
----	----	----	----	----	----	----	----	---	---	-----	---	-----	---	-----	------	---	------	-----	------	---	------	-----	------	---	------	-----	------	---	------	-----	------	---	------	-----	------	------

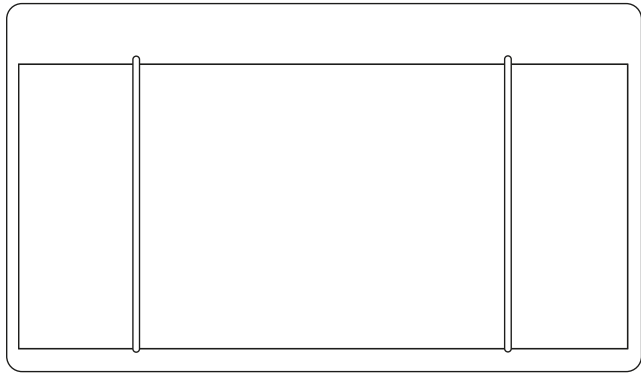

476 mm (18.7 in)

578 mm (22.8 in)

Trial Lens Set Comfort S Coated

The classic

The following modules are already included in the Trial Lens Set Comfort S coated

Description	Article no.	Quantity
tls trial lenses +/- sph (basic), coated	10033240	80
tls trial lenses +/- sph (extension), coated	10033242	56
tls trial lenses - cyl (basic), coated <i>axis position 90°</i>	10033248	20
tls trial lenses - cyl (extension), coated <i>axis position 90°</i>	10033252	16
tls trial lenses special lenses (basic) <i>2 x cross hair, 1 x occluder, 1 x pinhole disc 1.25</i>	10033256	4
tls trial lenses special lenses (extension) <i>coloured lens red, coloured lens green, frosted lens, plano lens, stenopaic slit 1.25 mm, Maddox cylinder red</i>	10033257	6
Trial Lens Set Comfort S, coated	10033419	182

We recommend the following optional module for expansion

Description	Article no.	Quantity
tls trial lenses +/- sph, coated, <i>each 2 x 16/18/20 D</i>	10033246	12

Trial Lens Trays Comfort S and Premium S allow you to combine the following optional modules:

	Options 1-7 (option 1 shown on the page before)	Article no.	Quantity
1	tls trial lenses + cyl (basic), coated, <i>axis position 90°</i>	10033250	20
	tls trial lenses + cyl (extension), coated, <i>axis position 90°</i>	10033254	16
	tls trial lenses prism. (basic), coated, <i>base 0°</i>	10033260	5
	tls trial lenses prism. (extension), coated, <i>base 0°</i>	10033262	5
2	tls trial lenses near-vision addition trial lenses, coated	10033258	22
	tls trial lenses prism. BO (basic/comfort), coated, <i>45°/135°</i>	10033263	24
3	tls trial lenses near-vision addition trial lenses, coated	10033258	22
	tls trial lenses prism. BI (basic/comfort), coated, <i>225°/315°</i>	10033264	24
4	tls trial lenses prism. BO (basic/comfort), coated, <i>45°/135°</i>	10033263	24
	tls trial lenses prism. BI (basic/comfort), coated, <i>225°/315°</i>	10033264	24
5	tls trial lenses near-vision addition trial lenses, coated	10033258	22
	tls trial lenses prism. (basic), coated, <i>base 0°</i>	10033260	5
	tls trial lenses prism. (extension), coated, <i>base 0°</i>	10033262	5
6	tls trial lenses prism. (basic/comfort), coated, <i>0°/180°</i>	10033268	24
	tls trial lenses prism. (extension/comfort), coated, <i>0°/180°</i>	10033269	12
7	tls trial lenses near-vision addition trial lenses, coated	10033258	22
	tls trial lenses prism. (basic/comfort), coated, <i>0°/180°</i>	10033268	24



The illustration above gives a schematic overview of the trial lens case (article no. 10035809) with the Trial Lens Set Comfort S, coated and the optional modules.

20	sph	20
18	18	18
16	16	16
14	14	14
13	13	13
12	12	12
11	11	11
10	10	10
9	9	9
8	8	8
7.5	7.5	7.5
7	7	7
6.5	6.5	6.5
6	6	6
5.5	5.5	5.5
5.25	5.25	5.25
5	5	5
4.75	4.75	4.75
4.5	4.5	4.5
4.25	4.25	4.25
4	4	4
3.75	3.75	3.75
3.5	3.5	3.5
3.25	3.25	3.25
3	3	3
2.75	2.75	2.75
2.5	2.5	2.5
2.25	2.25	2.25
2	2	2
1.75	1.75	1.75
1.5	1.5	1.5
1.25	1.25	1.25
1	1	1
0.75	0.75	0.75
0.5	0.5	0.5
0.25	0.25	0.25
0.12	0.12	0.12

6	cyl	6
5	5	5
4.5	4.5	4.5
4	4	4
3.5	3.5	3.5
3	3	3
2.75	2.75	2.75
2.5	2.5	2.5
2.25	2.25	2.25
2	2	2
1.75	1.75	1.75
1.5	1.5	1.5
1.25	1.25	1.25
1	1	1
0.75	0.75	0.75
0.5	0.5	0.5
0.25	0.25	0.25
0.12	0.12	0.12
MR	FL	
⊕	⊕	
⊖	⊖	
RF	GF	

3	sph	3
2.75	2.75	2.75
2.5	2.5	2.5
2.25	2.25	2.25
2	2	2
1.75	1.75	1.75
1.5	1.5	1.5
1.25	1.25	1.25
1	1	1
0.75	0.75	0.75
0.5	0.5	0.5
Δ	6	
5	5	5
4	4	4
3.5	3.5	3.5
3	3	3
2.5	2.5	2.5
2	2	2
1.5	1.5	1.5
1	1	1
0.75	0.75	0.75
0.5	0.5	0.5
0.25	0.25	0.25

20	sph	20
18	18	18
16	16	16
14	14	14
13	13	13
12	12	12
11	11	11
10	10	10
9	9	9
8	8	8
7.5	7.5	7.5
7	7	7
6.5	6.5	6.5
6	6	6
5.5	5.5	5.5
5.25	5.25	5.25
5	5	5
4.75	4.75	4.75
4.5	4.5	4.5
4.25	4.25	4.25
4	4	4
3.75	3.75	3.75
3.5	3.5	3.5
3.25	3.25	3.25
3	3	3
2.75	2.75	2.75
2.5	2.5	2.5
2.25	2.25	2.25
2	2	2
1.75	1.75	1.75
1.5	1.5	1.5
1.25	1.25	1.25
1	1	1
0.75	0.75	0.75
0.5	0.5	0.5
0.25	0.25	0.25
0.12	0.12	0.12

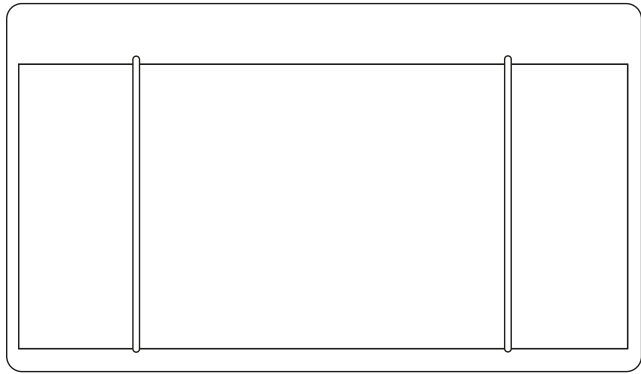

476 mm (18.7 in)

578 mm (22.8 in)

Trial Lens Set Premium S

The all-rounder

The following modules are already included in the Trial Lens Set Premium S coated

Description	Article no.	Quantity
tls trial lenses +/- sph (basic), coated	10033240	80
tls trial lenses +/- sph (extension), coated	10033242	56
tls trial lenses - cyl (basic), coated <i>axis position 90°</i>	10033248	20
tls trial lenses - cyl (extension), coated <i>axis position 90°</i>	10033252	16
tls trial lenses prism. BO (basic/comfort), coated <i>base 45°/135°</i>	10033263	24
tls trial lenses special lenses (basic) <i>2 x cross hair, 1 x occluder, 1 x pinhole disc 1.25</i>	10033256	4
tls trial lenses special lenses (extension) <i>coloured lens red, coloured lens green, frosted lens, plano lens, stenopaic slit 1.25 mm, Maddox cylinder red</i>	10033257	6
tls trial lenses near-vision addition trial lenses, coated	10033258	22
Trial Lens Set Premium S, coated	10033420	228



We recommend the following optional module for expansion

Description	Article no.	Quantity
tls trial lenses +/- sph, coated <i>each 2 x 16/18/20 D</i>	10033246	12

Everything perfectly organized!

You will find the ideal storage solution for your new trial lenses on **page 20**.

The illustration above gives a schematic overview of the trial lens case (article no. 10035809) with the Trial Lens Set Premium S, coated and the optional module.

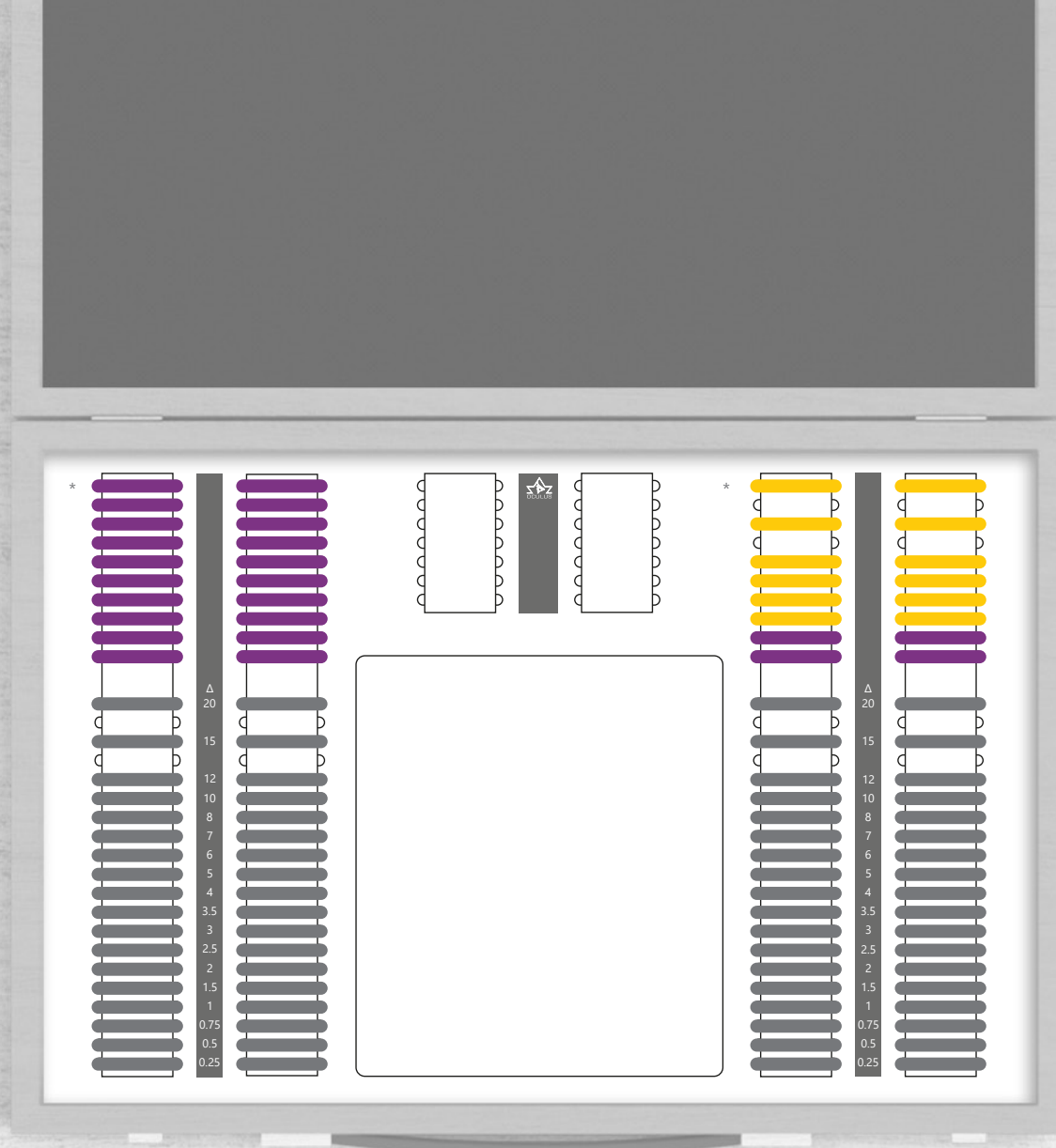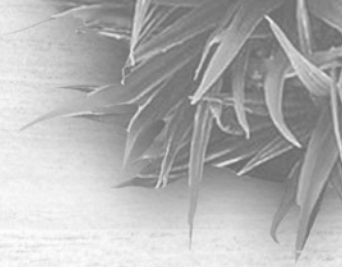

378 mm (14.9 in)

255 mm (10.0 in)



Trial Lens Set Prism S

The optimizer

The following modules are already included in the Trial Lens Set Prism S coated

Description	Article no.	Quantity
tls trial lenses prism. BO (basic/comfort), coated <i>base 45°/135°</i>	10033263	24
tls trial lenses prism. BI (basic/comfort), coated <i>base 225°/315°</i>	10033264	24
tls trial lenses prism. BO (extension/comfort), coated <i>base 45°/135°</i>	10033266	12
tls trial lenses prism. BI (extension/comfort), coated <i>base 225°/315°</i>	10033267	12
Trial Lens Set Prism S, coated	10033421	72

We recommend the following optional modules for expansion

Description	Article no.	Quantity
tls trial lenses prism. (basic/comfort), coated* <i>base 0°/180°</i>	10033268	24
tls trial lenses prism. (extension/comfort), coated* <i>base 0°/180°</i>	10033269	12

* Please note that the prism trial lens case has no diopters marking for these optional modules.

Everything perfectly organized!

You will find the ideal storage solution for your new trial lenses on **page 20**.

The illustration above gives a schematic overview of the trial lens case (article no. 10035806) with the Trial Lens Set Prism S, coated and the optional modules.

Overview of all Modules

Would you like to customize or supplement your trial lens case? Your distributor is ready to take your orders.

- 1 trial lens
- 2 trial lenses
- /•• 4 trial lenses
- BO base outside
- BI base inside

t/s trial lenses spherical [D]	+/- 20,00			••/••
	+/- 18,00			••/••
	+/- 16,00			••/••
	+/- 14,00		••/••	
	+/- 13,00		••/••	
	+/- 12,00		••/••	
	+/- 11,00		••/••	
	+/- 10,00		••/••	
	+/- 9,00		••/••	
	+/- 8,00	••/••		
	+/- 7,50		••/••	
	+/- 7,00		••/••	
	+/- 6,50		••/••	
	+/- 6,00	••/••		
	+/- 5,50		••/••	
	+/- 5,25		••/••	
	+/- 5,00	••/••		
	+/- 4,75		••/••	
	+/- 4,50	••/••		
	+/- 4,25		••/••	
	+/- 4,00	••/••		
	+/- 3,75	••/••		
	+/- 3,50	••/••		
	+/- 3,25	••/••		
	+/- 3,00	••/••		
+/- 2,75	••/••			
+/- 2,50	••/••			
+/- 2,25	••/••			
+/- 2,00	••/••			
+/- 1,75	••/••			
+/- 1,50	••/••			
+/- 1,25	••/••			
+/- 1,00	••/••			
+/- 0,75	••/••			
+/- 0,50	••/••			
+/- 0,25	••/••			
+/- 0,12		••/••		
coated	basic (10033240)	extension (10033242)	16/18/20 (10033246)	
uncoated	basic (10033238)	extension (10033241)	16/18/20 (10033243)	
quantity	80	56	12	



axis position 90° axis position 90°

t/s trial lenses - cylindrical [D]	- 6,00		••	
	- 5,00		••	
	- 4,50		••	
	- 4,00		••	
	- 3,50		••	
	- 3,00	••		
	- 2,75		••	
	- 2,50	••		
	- 2,25		••	
	- 2,00	••		
	- 1,75	••		
	- 1,50	••		
	- 1,25	••		
	- 1,00	••		
	- 0,75	••		
	- 0,50	••		
	- 0,25	••		
	- 0,12		••	
	coated	basic (10033248)	extension (10033252)	
	uncoated	basic (10033247)	extension (10033251)	
quantity	20	16		



axis position 90° axis position 90°

t/s trial lenses + cylindrical [D]	+ 6,00		••
	+ 5,00		••
	+ 4,50		••
	+ 4,00		••
	+ 3,50		••
	+ 3,00	••	
	+ 2,75		••
	+ 2,50	••	
	+ 2,25		••
	+ 2,00	••	
	+ 1,75	••	
	+ 1,50	••	
	+ 1,25	••	
	+ 1,00	••	
	+ 0,75	••	
	+ 0,50	••	
+ 0,25	••		
+ 0,12		••	
coated	basic (10033250)	extension (10033254)	
uncoated	basic (10033249)	extension (10033253)	
quantity	20	16	



t/s trial lenses prismatic [pdpt]	20,00			••		••		••
	15,00			••		••		••
	12,00		•	••		••		••
	10,00		•	••		••		••
	8,00		•	••		••		••
	7,00			••		••		••
	6,00		•	••		••		••
	5,00	•		••		••		••
	4,00	•		••		••		••
	3,50			••		••		••
	3,00	•		••		••		••
	2,50			••		••		••
	2,00	•		••		••		••
	1,50			••		••		••
	1,00	•		••		••		••
0,75			••		••		••	
0,50		•	••		••		••	
0,25			••		••		••	
coated	basic (10033260)	extension (10033262)	BO basic/comfort (10033263)	BO extension/comfort (10033266)	BI basic/comfort (10033264)	BI extension/comfort (10033267)	basic/comfort (10033268)	extension/comfort (10033269)
uncoated	basic (10033259)	extension (10033261)						
quantity	5	5	24	12	24	12	24	12

t/s trial lenses special lenses	Maddox cylinder red		•
	stenopaic slit 1.25		•
	plano lens, coated		•
	frosted lens		•
	coloured lens green		•
	coloured lens red		•
	pinhole disc 1.25	•	
	occluder	•	
	cross hair	••	
		basic (10033256)	extension (10033257)
quantity	4	6	

t/s trial lenses near-vision addition trial lenses [D]	+ 3,00	••
	+ 2,75	••
	+ 2,50	••
	+ 2,25	••
	+ 2,00	••
	+ 1,75	••
	+ 1,50	••
	+ 1,25	••
	+ 1,00	••
	+ 0,75	••
	+ 0,50	••
	coated	near-vision addition trial lenses (10033258)
	quantity	22

Storage – Trial Lens Cases and Lens Trays

Find the ideal solution for storing your trial lenses.



Trial Lens Case Basic S
Trial Lens Case Prism S

Beechwood case with handle, empty, with lens tray for 134 tls trial lenses

Dimensions (W x D x H):

378 x 255 x 55 mm (14.9 x 10.0 x 2.2 in)

- Article no. 10035803 (Trial Lens Case Basic S)
- Article no. 10035806 (Trial Lens Case Prism S)



Trial Lens Case Comfort S and Premium S / wall or table version

Trial Lens Case, empty, with lens tray for 244 tls trial lenses

Dimensions (W x D x H):

578 x 153 x 476 mm (22.8 x 6.0 x 18.7 in), wall version

- Article no. 10035812²⁾

Dimensions (W x D x H):

578 x 476 x 153 mm (22.8 x 18.7 x 6.0 in), table version

- Article no. 10035809²⁾



Lens tray for Trial Lens Set Basic S
Lens tray for Trial Lens Set Prism S

Lens tray, incl. diopters marking for 134 tls trial lenses

Dimensions (W x D x H):

357 x 234 x 23 mm (14.1 x 9.2 x 0.9 in)

- Article no. 10035813 (Lens tray for Trial Lens Set Basic S)
- Article no. 10035816 (Lens tray for Trial Lens Set Prism S)



Lens tray for Trial Lens Set Comfort S and Premium S

Lens tray, incl. diopters marking for 244 tls trial lenses

Dimensions (W x D x H):

545 x 401 x 36 mm (21.5 x 15.8 x 1.4 in)

- Article no. 10035814⁴⁾

OCULUS UB 5

The Universal Trial Frame for small heads and children

Whether for performing refractions on persons with a small or narrow head, myopia management or paediatric optometry – the UB 5 is flexible in use, making it the ideal solution for special requirements.

Highlights

- Five slots for each side: three in front, two behind
 - Gradation: 5°
 - PD adjustment R/L: 40 - 72 mm
 - Weight: 48 g
 - Parallax-free reading of the vertex distance (VD)
 - Flexible ear pieces
 - Anatomically designed nose pad for more comfort
 - Height- and depth-adjustable nose pad
 - Delivery content includes two additional nose pads
 - Thumb wheel for inclination adjustment
 - Adjustable arm length
 - Made of light metal/plastic
 - Delivery content includes five additional coloured knobs
- Article no. 42550




Info

The OCULUS UB 5 is the trial frame for trial lenses with a rim diameter of **28 mm** for the **modular concept tls**.

Trial Lens Cases

Technical Data

Trial Lens Case Basic S, empty (article no. 10035803)	
Dimensions (W x D x H)	378 x 255 x 55 mm (14.9 x 10.0 x 2.2 in)
Weight	approx. 2 kg (4.4 lbs)
Trial Lens Case Comfort S and Premium S, table version, empty (article no. 10035809)	
Dimensions (W x D x H)	578 x 476 x 153 mm (22.8 x 18.7 x 6.0 in)
Weight	approx. 6.8 kg (15 lbs)
Trial Lens Case Comfort S and Premium S, wall version, empty (article no. 10035812)	
Dimensions (W x D x H)	578 x 153 x 476 mm (22.8 x 6.0 x 18.7 in)
Weight	approx. 7 kg (15.4 lbs)
Trial Lens Case Prism S, empty (article no. 10035806)	
Dimensions (W x D x H)	378 x 255 x 55 mm (14.9 x 10.0 x 2.2 in)
Weight	approx. 2 kg (4.4 lbs)

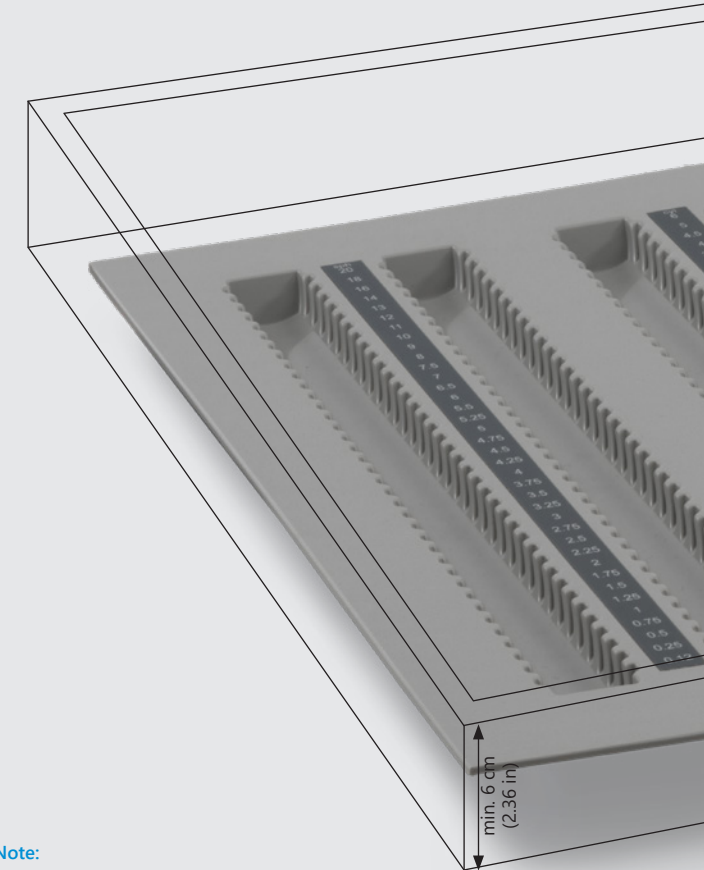

Note:

A drawer for keeping the lens tray requires a minimum height of 60 mm (2.36 in).

WWW.OCULUS.DE



The OCULUS QM system is certified in accordance with ISO 13485 (MDSAP) and (EU) 2017/745 (MDR)

OCULUS Optikgeräte GmbH
 Postfach • 35549 Wetzlar • GERMANY
 Tel. +49 641 2005-0 • Fax +49 641 2005-295
 Email: export@oculus.de • www.oculus.de

Find your local OCULUS representative on our website.

The availability of products and features may vary by country. OCULUS reserves the right to change product specifications and design. All information is valid at the time of printing (02/25).

05/0225/EN